


PCT

WORLD INTELLECTUAL PROPERTY ORGANIZATION  
International Bureau

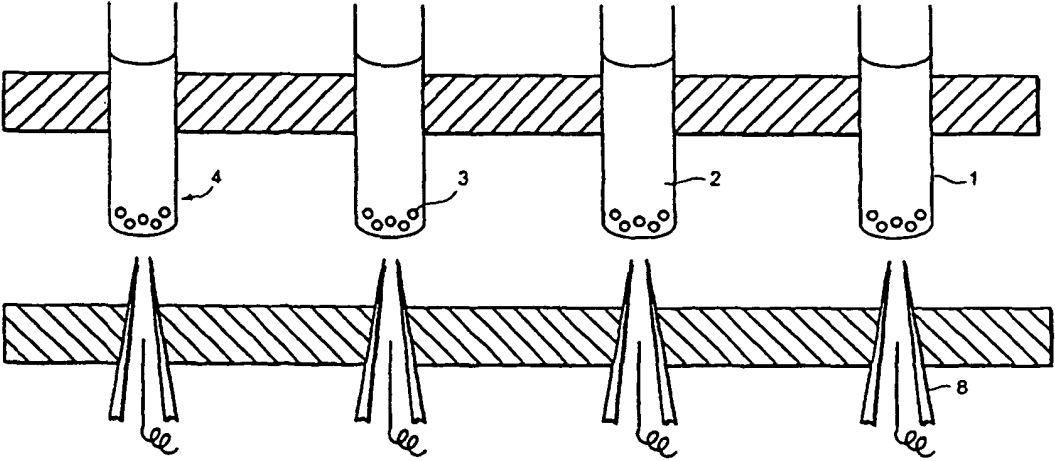

INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(51) International Patent Classification <sup>7</sup> : <b>G01N 33/487, 35/00, 35/10 // C12M 1/34, C12N 13/00</b>		A1	(11) International Publication Number: <b>WO 00/34776</b>
			(43) International Publication Date: <b>15 June 2000 (15.06.00)</b>
(21) International Application Number: <b>PCT/GB99/04073</b>			(81) Designated States: AE, AL, AM, AT, AU, AZ, BA, BB, BG, BR, BY, CA, CH, CN, CR, CU, CZ, DE, DK, DM, EE, ES, FI, GB, GD, GE, GH, GM, HR, HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MA, MD, MG, MK, MN, MW, MX, NO, NZ, PL, PT, RO, RU, SD, SE, SG, SI, SK, SL, TJ, TM, TR, TT, TZ, UA, UG, US, UZ, VN, YU, ZA, ZW, ARIPO patent (GH, GM, KE, LS, MW, SD, SL, SZ, TZ, UG, ZW), Eurasian patent (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM), European patent (AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE), OAPI patent (BF, BJ, CF, CG, CI, CM, GA, GN, GW, ML, MR, NE, SN, TD, TG).
(22) International Filing Date: <b>6 December 1999 (06.12.99)</b>			
(30) Priority Data: 9826742.0                      5 December 1998 (05.12.98)                      GB 9906053.5                      17 March 1999 (17.03.99)                      GB 9905998.2                      17 March 1999 (17.03.99)                      GB			
(71) Applicant (for all designated States except US): <b>CENES LIMITED [GB/GB]; Compass House, Vision Park, Chivers Way, Histon, Cambridge CB4 9ZR (GB).</b>			
(72) Inventors; and (75) Inventors/Applicants (for US only): <b>BYRNE, Nicholas, Gerard [GB/GB]; 31 Chivers Road, Haverhill, Suffolk CB9 9DS (GB). OWEN, David, Geraint [GB/GB]; 2 The Terrace, The Street, Kent DA12 3DF (GB).</b>			
(74) Agent: <b>DAVIES, Jonathan, Mark; Reddie &amp; Grose, 16 Theobalds Road, London WC1X 8PL (GB).</b>			<b>Published</b> <i>With international search report.</i>
(54) Title: <b>INTERFACE PATCH CLAMPING</b>			
			
(57) Abstract			
<p>The invention provides a novel development of the conventional patch clamp technique for measurement of whole cell electrical activity. The invention provides for one or more cell or cells to be suspended in a liquid medium at a liquid/air interface (by virtue of the effect of surface tension at the interface) whereby the cell or cells are accessible at the interface to a microstructure electrode (such as a pipette tip) to which a cell can attach to form an electrical seal, for the purpose of whole cell voltage clamp recording. According to the invention the electrode can be caused to form a high resistance electrical seal with a cell suspended in the liquid at the liquid/air interface without the need to press the cell against a solid support surface. The invention also provides apparatus for carrying out the interface patch clamp technique and control logic for operating a computer to carry out the interface patch clamp technique.</p>			

**FOR THE PURPOSES OF INFORMATION ONLY**

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AL	Albania	ES	Spain	LS	Lesotho	SI	Slovenia
AM	Armenia	FI	Finland	LT	Lithuania	SK	Slovakia
AT	Austria	FR	France	LU	Luxembourg	SN	Senegal
AU	Australia	GA	Gabon	LV	Larvia	SZ	Swaziland
AZ	Azerbaijan	GB	United Kingdom	MC	Monaco	TD	Chad
BA	Bosnia and Herzegovina	GE	Georgia	MD	Republic of Moldova	TG	Togo
BB	Barbados	GH	Ghana	MG	Madagascar	TJ	Tajikistan
BE	Belgium	GN	Guinea	MK	The former Yugoslav Republic of Macedonia	TM	Turkmenistan
BF	Burkina Faso	GR	Greece	ML	Mali	TR	Turkey
BG	Bulgaria	HU	Hungary	MN	Mongolia	TT	Trinidad and Tobago
BJ	Benin	IE	Ireland	MR	Mauritania	UA	Ukraine
BR	Brazil	IL	Israel	MW	Malawi	UG	Uganda
BY	Belarus	IS	Iceland	MX	Mexico	US	United States of America
CA	Canada	IT	Italy	NE	Niger	UZ	Uzbekistan
CF	Central African Republic	JP	Japan	NL	Netherlands	VN	Viet Nam
CG	Congo	KE	Kenya	NO	Norway	YU	Yugoslavia
CH	Switzerland	KG	Kyrgyzstan	NZ	New Zealand	ZW	Zimbabwe
CI	Côte d'Ivoire	KP	Democratic People's Republic of Korea	PL	Poland		
CM	Cameroon	KR	Republic of Korea	PT	Portugal		
CN	China	KZ	Kazakstan	RO	Romania		
CU	Cuba	LC	Saint Lucia	RU	Russian Federation		
CZ	Czech Republic	LI	Liechtenstein	SD	Sudan		
DE	Germany	LK	Sri Lanka	SE	Sweden		
DK	Denmark	LR	Liberia	SG	Singapore		

- 1 -

## INTERFACE PATCH CLAMPING

## Introduction

The present invention provides a novel development of the conventional patch clamp technique. This novel technique  
5 is referred to as the interface patch clamp method.

Voltage gated ion channels are potential targets for a considerable range of novel treatments in a variety of disease states. The development of the patch clamp technique has provided a powerful method for the study of  
10 ion channel function and pharmacology in whole cells. However, while the patch clamp technique provides a definitive method for the investigation and screening of drugs with potential activity on voltage gated ion channels, the technique is currently highly dependent on  
15 the skill of the operator and tends to be very slow for drug screening. The present invention provides a method for increasing the rate at which compounds may be screened for ion channel blocking/agonist activity using the patch clamp technique. The method can retain the essential  
20 features of the conventional patch clamp recording system while facilitating automation of the major time-consuming components of the technique.

## Background: Conventional Patch Clamp

The success of the patch clamp technique is derived from  
25 the ability to form "tight" (i.e. high resistance: Giga Ohm) electrical seals between an area of the cell membrane (the Patch) and the tip of a pipette. The patch clamp pipette is usually made from glass. The formation of the G-seal is dependent on the profile of the top of the  
30 pipette, and is enhanced by the application of suction to

- 2 -

the interior of the pipette. The requirements for the formation of the G-seals are well established and the process is usually monitored electrically by display of the current pulse recorded in response to a small voltage step applied throughout seal formation. After formation of a G-seal, the area of membrane under the pipette may be disrupted to obtain whole cell voltage clamp recording mode.

The sequence of events leading to successful G-seal formation and whole cell recording mode using pre-formed patch pipettes is as follows:

1. Selection of a suitable cell.
2. The patch pipette is positioned approximately 50 microns above the cell.
3. The pipette is lowered until the cell surface is deformed by the pipette tip.
4. Negative pressure is applied to the interior of the pipette until a G-seal is formed between the pipette tip and the cell membrane.
5. Whole cell recording mode is established by the application of further negative pressure which disrupts the cell membrane in the area under the pipette tip.

Steps two and three are slow and require considerable manual dexterity and a high level of operator skill. Visualisation of the cells and the patch pipette requires the use of a high quality microscope and, in order to position the pipette, a high quality three axis

- 3 -

micromanipulator with sub-micron resolution in each axis is required.

#### Summary of the Invention

In its broadest terms the invention provides for one or  
5 more cell or cells to be suspended in a liquid medium at a liquid/air interface (by virtue of the effect of surface tension at the interface) whereby the cell or cells are accessible at the interface to a microstructure electrode (such as a pipette tip) to which a cell can attach to form  
10 an electrical seal, for the purpose of whole cell voltage clamp recording. According to the invention the electrode can be caused to form a high resistance electrical seal with a cell suspended in the liquid at the liquid/air interface without the need to press the cell against a  
15 solid support surface.

Any body of liquid or column of liquid, which gives rise to a situation in which a cell or cells are located in the liquid at a liquid/air interface, can be used in the invention. For examples cells may be suspended in a  
20 column of liquid held by surface tension in a capillary tube. Alternatively cells may be suspended in a droplet of liquid, which droplet may itself be suspended from or supported by a support.

It will readily be appreciate that the interface patch  
25 clamp technique can be operated in "single cell mode", or could be multiplexed to operate on a matrix of cells with multiple electrodes.

According to one aspect of the invention, interface patching can utilise a patch pipette of conventional type.  
30 Cells are supported on a liquid/air interface at one end

- 4 -

of a capillary tube (e.g. made of glass, polyethylene or other suitable material). The axis of the patch pipette is in line with the axis of the tube so that the pipette tip can be manipulated into the opening of the tube where the cells are supported at the air/liquid interface. The capillary tube or the patch pipette can be mounted onto a single axis manipulator. Only one manipulator is required and this may be used to move either the patch pipette or the capillary tube. Whole cell recording mode is established as follows:

6. A layer of cells is established at the interface between the extracellular physiological solution (the liquid in which the cells are suspended) and air by dipping the capillary tube into a suspension of cells. The density of cells in the suspension must be sufficient to provide a sufficient number of cells to form a layer of cells at the interface.
7. Electrical contact with the extracellular solution is established via a non-polarizable electrode (e.g. an Ag/AgCl wire) and the tube is mounted either to a fixed clamp or single axis manipulator.
8. A patch pipette is provided which can be filled with electrolyte solution.
9. The patch pipette is mounted concentrically with the capillary tube either via a single axis manipulator or fixed clamp (if the capillary tube is to be moved). The pipette filling solution is connected via the non-polarizable electrode to the headstage of a conventional patch clamp amplifier. The pipette holder allows suction to be applied to the pipette interior.

- 5 -

10. Cell attached patch mode of recording is established by bringing the pipette tip in contact with the interface by moving the pipette and the capillary tube respectively together along the single mounting axis (e.g. either by moving the pipette towards the tube and interface or vice versa). On entry into the interface the movement of the pipette and capillary tube together is stopped and the pipette current is offset to zero on the patch clamp amplifier. The resistance of the pipette increases when the pipette contacts one of the cells at the air/liquid interface. Suction is then applied to the interior of the pipette and the pipette and capillary tube are moved closer together until the pipette tip is located inside the capillary tube.

Initial seal formation between the pipette tip and the cell may also be assisted by the application of gentle suction during entry of the pipette into the interface.

- 20 A G-seal is formed between the patch pipette tip and the cell membrane by the application of further suction to the interior of the pipette and monitoring the pipette resistance.

11. Following the formation of cell attached patch mode, the suction is released, pipette current is offset to zero and a holding voltage applied to the pipette (e.g. -60mV).

12. A whole cell recording is obtained by the application of further suction to the pipette interior until the whole cell recording mode is established in conventional manner.

- 6 -

According to this invention it is preferred that the capillary tube should be mounted in an upright orientation (i.e. essentially vertically) with the air/liquid interface at the downward end of the tube.

5 This has the advantage that suspended cells will tend to "sediment" naturally to the downward end of the tube and be collected there in a layer. The layer will preferably be several cells deep and loosely packed. Thus according to the invention the pipette tip may be moved upwardly  
10 relative to the air/liquid interface at the tube end (either by moving the pipette or the tube along the single axis) so as to come into contact with a cell in the layer at the interface. The relative density or concentration of cells at the interface compared to the density in the  
15 bulk of the liquid in the tube ensures a high probability that a cell can be collected on the tip without the need for visualisation of the operation and without the need for multidirectional manipulation of the tip/cell positional relationship. Surprisingly it has been found  
20 that G-seal formation between the cell and the pipette can occur without pressing the cell against a solid substrate.

Where the arrangement is intended to operate with the pipette in an upright orientation (i.e. essentially vertically) with the tip uppermost and pointing upwardly,  
25 the pipette should be constructed so as to prevent the filling electrolyte solution flowing out and being lost. This may be achieved for example by use of a custom-made mounting assembly and/or by shaping the pipette body to prevent loss of filling solution (e.g by bending the  
30 pipette shaft into a U- or J- shape).



- 7 -

The invention also provides methods and apparatus employing control logic to allow automation of a patch clamp system employing the Interface Patch Clamp technique described herein. The logic described will control one or more electromechanical micromanipulators/translators holding one or more patch clamp pipettes and/or capillary tubes in order to patch clamp cells and apply drugs/compounds in order to screen for activity on membrane ion channels. A major advantage of the logic described is that automation is achieved in this system by the use of feedback from signals from the patch clamp amplifier and no image recognition software is required.

The invention is illustrated by way of example in the accompanying figures in which:

- Figure 1a shows a capillary tube containing a suspension of cells; and  
Figure 1b shows the cells having formed a layer at the air/liquid interface at one end of the capillary tube;  
Figure 2 shows a general arrangement of the interface patch clamp recording equipment with moveable capillary tube;  
Figure 2a shows an Apparatus for Interface Patch Clamping with drug/compound application;  
Figure 3 shows the cell attached to the patch pipette ready for recording mode.  
Figure 4 shows drug/compound addition during interface patch clamp recording: start position;  
Figure 5 shows drug/compound addition during interface patch clamp recording: extracellular solution added to dish and dish moved down;  
Figure 6 shows drug/compound addition during interface patch clamp recording: solution in dish brought into contact with interface region;

- 8 -

Figure 7 shows drug/compound addition during interface patch clamp recording: capillary raised above surface of solution in dish;

5 Figure 8 is a flow diagram of control logic embodying a further aspect of the invention;

Figure 9 is a flow diagram of an example of the G-seal formation steps of figure 8;

Figures 10 to 16 are flowcharts of a third embodiment of the control logic used in the invention, wherein

10 Figure 10 is a flow diagram of control logic embodying a further aspect of the invention;

Figure 10a is a flow diagram expanding the reset routine of figure 10;

15 Figure 10b is a flow diagram expanding the safe quit routine of figure 10;

Figure 11 is a flow diagram expanding the interface or cell locator routine of figure 10;

Figure 11a is a flow diagram expanding the interface hit routine of figure 11;

20 Figure 12 is a flow diagram expanding the Giga seal test routine of figure 10;

Figure 13 is a flow diagram expanding the whole cell detection routine of figure 10;

25 Figure 14 is a flow diagram expanding the qualification routine of figure 10;

Figure 15 is a flow diagram expanding the quality monitor routine of figure 14;

Figure 16 is a flow diagram expanding the experiment routine of figure 10;

30 Figure 16a is a flow diagram expanding the run trace routine of figure 16;

Figure 16b is a flow diagram expanding the stabilise routine of figure 16;

35 Figures 17 to 20 show recordings of current (1) and voltage (2) obtained from an automated patch clamp system

- 9 -

(AutoPatch) under software control using the Interface Patch Clamp technique. (Recording from an MK1 cell); Figure 21 shows the effect of the potassium channel blocking drug tetraethylammonium (TEA) on the potassium current recorded from MK1 cell in whole cell recording mode obtained using the Interface Patch Clamp technique; Figure 22 shows cells suspended in a droplet of liquid covering a hole through a support surface and providing an air/liquid interface according to the invention; Figure 23 shows a further embodiment of droplets containing cell suspensions forming air/liquid interfaces according to the invention; Figure 24 shows a multiplexed interface patch clamping array; Figure 25 shows an alternative form of electrode microstructure according to the invention; and Figures 26a and b show a multiwell arrangement for drug application to patch-clamped cells.

Referring to Figure 1a; a capillary tube (1) of appropriate size can pickup and hold a liquid sample (2) containing cells (3) in suspension. The sample can be picked up simply by dipping the tube end into a suitable bulk liquid reservoir. The liquid in the tube forms an air/liquid interface (4) at the tube end (5). The cells are initially distributed throughout the liquid relatively evenly.

Referring to Figure 1b; with the tube in an upright essentially vertical orientation, the cells tend to sediment and to pack loosely together at the lower end of the tube by the tube end to form a layer (6) several cells deep. It will be appreciated by those skilled in the art that the density and depth of the cell layer can be determined by such factors as the cell concentration in

- 10 -

the original suspension, the sedimentation time, the relative density of the cells and the liquid etc. It will also be appreciated that means could be devised to encourage or assist cells to migrate from the liquid  
5 towards the air/liquid interface rather than or as well as relying on gravitational sedimentation alone. The Figure also shows the top of a patch pipette 8 pointing upwardly towards the interface.

Referring to Figure 2; an arrangement is shown in which a  
10 single axis manipulator is used to move a capillary tube 1 held in a clamp (7) relative to a fixed patch pipette (8) held in a clamp (9). It will be apparent to those skilled in the art that this could be reversed so that the pipette is moved and the tube is fixed. The Figure shows the tube  
15 clamped in a linear bearing sliding block (10) attached to a motorised single axis manipulator (11). The manipulator should be controlled preferably by computer in order to allow the motion of the manipulator to be varied by feedback from the patch clamp amplifier. The patch  
20 pipette is provided with a connection (12) to a conventional headstage. The system is also provided with a source of variable suction under the control of the patch clamp amplifier/computer.

In Figure 2a an arrangement is shown in which additional  
25 electromechanical micromanipulators have been added. The micromanipulator labelled (13) is for moving the glass capillary under automated or manual control. A second micromanipulator (14) moves the dish for drug application up and down the glass capillary. A third micromanipulator  
30 (15) moves a modified pipette holder to provide electrical contact with the pipette and a means of applying suction to the interior of the pipette. Rotational bases (16 and 17) allow the pipette holder to be moved in and out of the

- 11 -

recording area and rotation of the pipette through 180 degrees for filling with pipette solution.

The Figure also shows additional features, namely; a pipette holder (18); a patch clamp headstage (19); and a  
5 dish holder (20).

A version of the apparatus is envisaged in which patch pipettes will be loaded and filled automatically under software control. It is envisaged also that the loading of capillary glass into the apparatus and the filling with  
10 cell suspension will also be automated.

Referring to Figure 3; a G-sealed cell 3 is shown held on the tip of the patch pipette 8 and positioned within the entrapped liquid volume in the tube.

Cell attached patch and whole cell (voltage clamp)  
15 recording may then be carried out.

The invention described herein has a number of significant features:

- Visualisation of the pipette and the cell is not  
20 required.
- Novel recording configuration that would not be considered as obvious.
- Surprisingly G-seal formation occurs without pressing the cell against a hard substrate.
- 25 • Cells form a layer at the solution-air interface.
- G-seal formation may be achieved using electronic feedback alone.
- There is no requirement for optical recognition/feedback.
- 30 • The system can be automated.

- 12 -

- Multiple recording capillaries and pipettes may be employed in order to allow recordings to be made simultaneously from many cells.

Exemplary methods of operation of the apparatus of the  
5 embodiment under software control to achieve various of  
these advantages are described below.

In order to use the invention for screening compound (e.g.  
for ion channel blocking/agonist activity) the compound of  
interest needs to be applied to the cell attached to the  
10 patch pipette. It will readily be appreciated that this  
could be achieved in different ways, for example by adding  
the compound to the extracellular liquid in the capillary  
tube either before or after G-seal formation. One  
additional advantage of the invention is that the liquid  
15 in the tube could be arranged in layers (e.g. containing  
different compounds or different concentrations of  
compounds) and the single axis manipulator could then be  
used to physically move and position a cell on a pipette  
tip into a chosen layer (e.g. by moving the G-sealed cell  
20 on the tip further up the tube away from the air/liquid  
interface at one tube end).

A further example of how the effects of compounds may be  
studied is illustrated in figures 4 to 7.

Figure 4 shows a capillary (1) containing the cell  
25 suspension (2) and patch pipette (8) in the recording  
position for whole cell recording from a cell at the  
pipette tip. In addition, the capillary tube has been  
inserted through a hole (21) made in a dish (22) (e.g.  
35mm plastic culture dish or similar). The dish is made  
30 of a material with hydrophobic properties and the hole

- 13 -

allows the dish to be raised and lowered along the axis of the capillary by means of a micromanipulator (14).

Figure 5 shows the dish after it has been filled with extracellular physiological solution (23), which may contain the drug to be studied, or the drug may be added at a later stage. Surprisingly, if the fluid level in the dish is low, leakage through the hole does not occur because the tendency to leak is counterbalanced by:

1. The surface tension of the water
2. The attraction of the water/solution to the glass capillary.

After adding the solution to the dish, it is lowered in the direction of the arrow.

Figure 6 shows the solution in the dish in contact with the end of the glass capillary and the patch pipette. The dish and the capillary are now raised simultaneously (arrows A and B) in order to position the pipette tip/cell within the layer of liquid in the dish. If drug is present in the dish at this point and the capillary and dish were moved upwards rapidly, this would constitute a rapid application system particularly useful for the study of agonist responses that desensitise.

Figure 7 shows the effect of raising the capillary so that it is not in contact with the liquid in the dish. The pipette tip/cell remains immersed in the external solution layer in the dish. The solution may be exchanged readily by perfusion of the dish and this allows multiple drug additions and dose response curves to be obtained while recording from the one cell.

- 14 -

## Control Logic for an Automated Patch Clamp System

### Introduction

The following describes three embodiments of the control logic required to allow automation of a patch clamp system employing the Interface Patch Clamp technique described herein. In each case, the logic described will control one or more electromechanical micromanipulators/translators in order to patch clamp cells and apply drugs/compounds in order to screen for activity on membrane ion channels. A major advantage of the logic described is that automation is achieved in this system by the use of feedback from signals from the patch clamp amplifier and no image recognition software is required.

### 15 Methods

*Inputs to the program in all the embodiments are required from the patch clamp amplifier as follows:*

Imon = current monitor output

Vhold = holding potential

20 *Inputs to the program derived from patch clamp amplifier output signals are required as follows:*

Inoise = base line current noise recorded from Imon

Rpip = pipette resistance

Rtot = Total resistance

25 Rs derived from Imon signal during voltage step

It is envisaged that these signals and evaluated values will be obtained from existing software (such as Heka



- 15 -

software) via a suitable software interface. These signals and evaluated values are further defined in the list of variables and parameters below.

*Inputs from manipulators/translators are required as follows, or from the following devices.*

- Patch module micromanipulator encoder
- Capillary clamp/loader encoder and empty signal
- Pipette automated clamp/loader encoder and empty signal
- Two axis translator encoders for cell dipper
- Drug application micromanipulator encoder
- Pipette holder micromanipulator encoder

*Control outputs from computer are required for the following devices.*

- Patch module manipulator
- Pipette automated holder
- Capillary loader/clamp
- Pipette loader/clamp

- 16 -

Two axis translator for cell loading system  
Pipette clamp  
Suction device  
Drug application manipulator  
Drug perfusion solenoid valve system

The software uses signals derived from the patch clamp amplifier in order to control a number of peripheral devices which carry out patch clamping using the Interface Patch Clamp technique. The devices controlled by the  
5 logic comprise a number of micromanipulators, a suction device for the patch pipette and a valve system for perfusion of a recording chamber, such as the dish described above.

A number of parameters are given pre-set values which can  
10 be changed by the operator to suit different experimental conditions.

Summary of the control logic for the automated Interface Patch Clamp - First Embodiment

Initial seal formation

15 The sequence of movements required for formation of a G-seal is unique for interface patch clamping and involves the control of at least one single axis manipulator (e.g. the patch module motor, although either the pipette or the capillary may be moved to achieve the necessary relative  
20 movement between them) with feedback from the patch clamp amplifier. In a first embodiment of the control logic, the pipette is initially spaced from the capillary, as illustrated in Figure 1b for example, and is moved towards the liquid/air interface at the capillary end until a  
25 change in the current monitor signal is recorded when the patch pipette enters the liquid/air interface and this

- 17 -

signal is used as the trigger to stop the micromanipulator. The pipette resistance may be derived from the output of the patch clamp amplifier and initial seal formation is monitored by recording the change in  
5 pipette resistance. If the resistance of the pipette does not increase beyond a pre-set value, the control logic infers that no G-seal has been formed and activates the patch module motor to move the liquid/air interface and the pipette apart until the resistance increases, which  
10 may occur when the pipette tip is withdrawn from the liquid or when a narrow neck of liquid is drawn out by surface tension between the pipette tip and the capillary end. When a resistance increase to a pre-set value is recorded suction is applied to the interior of the patch  
15 pipette and the patch pipette and the liquid/air interface are moved towards each other to a pre-set point, in a further attempt to form a G-seal with a cell.

#### Whole cell recording mode

After formation of the cell attached patch clamp recording  
20 mode, whole cell mode is obtained by the application of suction to the interior of the patch pipette while simultaneously monitoring the current ( $I_{\text{mon}}$ ) for capacitative transients. In the logic described, the formation of whole cell recording mode is detected by a  
25 threshold crossing method but it will be apparent to those skilled in the art that other methods may also be employed e.g. online FFT (Fast Fourier Transform), Template Matching etc. The control logic checks for incorrect  
30 detection of whole cell mode before activating the experimental protocol.

#### Cell quality test

- 18 -

This routine monitors the quality of voltage clamp by comparing the series recorded resistance with a value related to pipette resistance. It will be appreciated that this method may be further enhanced by relating the acceptable series resistance to the amplitude of current evoked by a voltage step. In addition, an additional loop may be added to include the possibility of recording with a maintained level of suction applied to the pipette in cells that exhibit continuously increasing values of series resistance. The quality of the cell is also monitored by the holding current which should not be more negative than a pre-set value. It will be appreciated that this method could be enhanced by relating the acceptable value for holding current to the amplitude of the current in response to a voltage step.

#### Drug/compound application

The initial phase of drug application is unique to the interface patch clamp technique and involves the control of two single axis micromanipulators. The movements required utilise the position of the patch module micromanipulator recorded on entry into the interface as a reference point. After the cell has been immersed into solution contained in a perfusion chamber, the control logic calls a routine to carry out perfusion of the chamber via the activation of solenoid flow control valves.

#### Control Logic in detail - (First Embodiment)

##### *Variables/parameters*

P = Pipette pressure relative to atmospheric pressure defined as 0  
d = Patch module motor position  
d0 = Patch module motor start position  
d1 = Patch module motor position following entry of pipette into interface  
d2 = Patch module motor position with pipette in recording position  
d3 = Patch module motor position for chamber perfusion  
dapp = Drug application module micromanipulator position  
dapp0 = Drug application module micromanipulator position 0  
dapp1 = drug application micromanipulator position increment 1 Pre-set increment  
dapp2 = drug application micromanipulator position increment 2 Pre-set increment  
Rs = Series resistance  
Cslow = Slow capacitance compensation  
Cfast = Fast capacitance compensation  
Inoise = base line noise  
Rpip = pipette resistance  
Rtot = Total resistance  
R1 Initial seal resistance (pre-set value)  
R2 Seal resistance required for progression to whole cell

- 19 -

Vhold = Pre set holding potential in mV  
 Imon = current monitor output  
 i = pre-set holding current  
 Ihold base line holding current for holding potential  
 dpip = pipette holder module motor position  
 dpip0 = pipette holder module motor position start  
 dpip1 = pipette holder module motor position pipette on  
 pclamp = position of pipette clamp/loader encoder  
 pclamp = 0 pipette not clamped (loading position)  
 pclamp = 1 pipette clamped (recording position)  
 rotclamp = position of rotory stage mounting for pipette loader  
 rotclamp = 0 recording position  
 rotclamp = 1 pipette filling position  
 pipfil = pipette filler position  
 pipfil = 0 pre-/post-fill position  
 pipfil = 1 fill position  
 pipsyringe = pipette filler driven syringe position  
 pipsyringe = dv driven syringe movement required to fill pipette (pre-set value)  
 cclamp = position of capillary clamp/loader encoder  
 cclamp = 0 capillary not clamped (loading position)  
 cclamp = 1 capillary clamped (recording position)  
 pload = pipette loader empty signal  
 pload = 0 pipettes in loader  
 pload = 1 pipette loader empty  
 cload = capillary loader empty signal  
 cload = 0 capillaries in loader  
 cload = 1 capillary loader empty  
 celldiph = horizontal translator for cell dip  
 celldiph = 0 cell storage encoder position (pre-set)  
 celldiph = 1 dip encoder position (pre-set)  
 celldipv = verticle translator for cell dip  
 celldipv = 0 pre-/post- dip encoder position (pre-set)  
 celldipv = 1 capillary dip encoder position (pre-set)  
 tdelay = variable delay between clamping capillary and starting to patch clamp  
 dt = time interval  
 dt1 = pre-set waiting time interval suction off (s)  
 dt2 = pre-set suction time interval (s)  
 dt3 = pre-set suction time interval (s)  
 dt4 = pre-set suction time interval (s)  
 dt5 = pre-set suction time interval (s)  
 x = suction increment factor  
 f = frequency of seal test pulse  
 detectmin = 0 -ve capacitance transient (3Inoise threshold) not detected  
 detectmin = 1 -ve capacitance transient (3Inoise threshold) detected  
 detectmax = 0 +ve capacitance transient (3Inoise threshold) not detected  
 detectmax = 1 +ve capacitance transient (3Inoise threshold) detected

- 20 -

I = Whole cell mode flag  
 I = 0 Not whole cell  
 I = 1 Whole cell mode established  
 singlemV = pre-set voltage test pulse  
 curr = current recorded between pre-set cursors during voltage step  
 testcurr = pre-set value for current required to start experimental protocol  
 Valve1 - 8 = solenoid valves controlling supply of solution to perfusion dish  
 tv = time interval for valve activation (pre-set)  
 drain valve = controls suction supply to perfusion dish

### Control Logic - Second Embodiment

Control logic according to a further embodiment is illustrated as a flow diagram in Figure 8. Exemplary logic steps within each of the functions shown in the flow diagram are set out below.

#### 00 Initialisation

d = d0  
 dpip = dpip0  
 pclamp = 0  
 pipfil = 0  
 rotclamp = 0  
 cclamp = 0  
 celldiph = 0  
 dapp = dapp0  
 Rtot  $\geq$  20M  
 Imon = Inoise  
 P = 0  
 If pload = 0 and cload = 0 then GOTO 01  
 If pload = 1 then report "Re-load pipette cassette" and GOTO 19  
 If cload = 1 then report "Re-load capillary cassette" and GOTO 19

#### 01 Autofeed

Move capillary clamp motor cclamp = 1  
 Move celldiph translator to celldiph = 1  
 Move celldipv translator to celldipv = 1  
 Move celldipv translator to celldipv = 0  
 Move celldiph translator to celldiph = 0  
 start timer  
 wait for variable delay = tdelay  
 GOTO 02

#### 02 Pipette load/fill

- 21 -

Move pipette clamp motor pclamp = 1  
Move rotation stage 180 degrees rotclamp = 1  
move pipfil until pipfil = 1  
move motor driven syringe drive until pipsyringe = dv  
move pipfil until pipfil = 0  
move rotation stage 180 degrees rotclamp = 0  
move dpip = dpip1  
GOTO 03

## 03 Junction null

Seal test signal on  
Compensate Cfast  
Move patch module down until Imon >/< Inoise  
5 Record patch module motor position d=d1  
Activate Junction null  
Measure Rtot  
Rtot=Rpip  
If Rpip <10M and >/=4M GOTO 04 else Report "Pipette  
10 resistance out of range" and GOTO 20

## 04 Formation of Gseal

Measure Rtot  
If Rtot=/> 2Rpip  
Suction on P=-pmmHg  
15 Move patch module down until d=d2  
(d2 =pre-set recording position)  
GOTO 05

## 04.1

If Rtot <2 Rpip  
20 Wait for time t1  
After time t1 move patch module upwards  
until Rtot>2Rpip then stop. Wait until

- 22 -

Rtot=2Rpip then move patch module down  
until d=d1

If d=d1 and Rtot=> 2Rpip then suction on  
P = -pmmHg

5 and move patch motor module down until  
d=d2, else GOTO 04.1 until maximum of 5  
iterations, then GOTO 20

05 Seal test loop

10           05.1           N=N+1



- 23 -

Measure  $R_{tot}$  with  $-pmmHg$  for time interval  $\Delta t_2$

If  $R_{tot} \geq R_1$  and  $\Delta t \leq \Delta t_2$

Then suction off  $P = P_0$  until  $R_{tot} \neq R_2$  or  $\Delta t = t_1$

If  $R_{tot} < R_1$  and  $\Delta t > t_1$

Then repeat 05.1 until  $N=5$  or  $R_{tot} \neq R_2$

If  $N = 5$  and  $R_{tot} < R_2$  Then suction on  $-xpmmHg$

Repeat 05.1 until  $N = 5$  or  $R_{tot} \neq R_2$

If  $N = 5$  and  $R_{tot} < R_2$  Then suction on  $-xpmmHg$

Repeat 05.1 until  $N = 5$  or  $R_{tot} \neq R_2$  or  $x = p_{max}$

If  $x = p_{max}$  and  $R_{tot} < R_2$  Then Report "Unable to obtain G seal" GOTO 20

If  $R_{tot} \neq R_2$  then GOTO 06 for Whole Cell Mode or  
GOTO 08 for Cell Attached Patch Mode

#### 06 Whole Cell - Threshold Method

Compensate  $C_{fast}$

$I = 0$

$h_p = v_{hold}$

06.1 Suction on  $-pmmHg$  until  $I=1$  or  $\Delta t = \Delta t_3$

06.2 Transient detection

$N = N + 1$

If  $detect_{max} = 1$  and  $detect_{min} = 1$

Then  $I = 1$

Repeat 06.2 until  $detect_{max} = 0$  and  $detect_{min} = 0$  or  $N = 10$

If  $N = 10$  Then GOTO 07

If  $detect_{max} = 0$  and  $detect_{min} = 0$

Then GOTO 06.1

#### 07 Cell quality test

Measure  $R_s$  and  $I_{hold}$

If  $R_s \neq 3R_{pip}$  Suction on  $-pmmHg$  and start timer

- 24 -

When  $R_s < 3R_{pip}$  suction off

If  $I_{hold} < i_{pA}$  then and  $R_s < 3R_{pip}$  then seal test signal off  
and GOTO 08

If  $I_{hold} > i_{pA}$  and time int = dt4 then GOTO 20

If  $R_s \neq 3R_{pip}$  and time int = dt5 then suction on  
-pmmHg

If  $R_s < 3R_{pip}$  and  $I_{hold} < i_{pA}$  then seal test signal off and  
GOTO 08

#### 08 Experimental protocol

Apply single voltage step to single mV  
Measure current amp during voltage steps  
If curr < testcurr then stop voltage steps  
Report "Control current out of range"  
GOTO 15

If curr  $\neq$  testcurr activate voltage step protocol GOTO 09

Measure  $R_s$   
If  $R_s \neq 3R_{pip}$  stop voltage protocol and GOTO 07  
Measure  $I_{hold}$  during interval between voltage steps

If  $I_{hold} > i_{pA}$  then stop voltage step/drug application  
protocol and GOTO 01

#### 09 Voltage step protocol

Uses program already available  
Program must call drug application sub-routine 10

#### 10 Drug/compound application

drain valve on  
Fill drug application chamber Valve1 on for time interval

- 25 -

Move drug application micromanipulator down to d1  
-dapp1 (NB -ve value represents downward movement)  
Move dapp = d1 + dapp2 and d= d2 + dapp2  
simultaneously  
Move d = d3 glass capillary moved up  
Call software controlling flow control valves

15 Reset Autopatch

GOTO 00.

19, 20 Steps 16 onwards are routines for setting  
up the apparatus controlled by the  
software/method and do not relate to the  
inventive operation of the apparatus and  
5 their design is within the normal  
competency of the skilled person. Steps  
19 and 20 relate respectively to reloading  
the pipette and capillary cassettes and to  
checking and/or resetting of the apparatus  
10 if operation is unsuccessful.

Figure 9 is a flow diagram illustrating steps 03 and 04,  
which relate to formation of the G-seal. These steps  
comprise the most important advantageous steps in this  
embodiment for controlling the method and apparatus  
15 described herein.

In step 03 (junction null), the pipette tip is initially  
spaced below the meniscus at the end of the capillary.  
The logic, or software, then controls the patch module  
motor to move the pipette tip towards the meniscus until  
20 contact is made, detected by electrical contact  
therebetween. The movement is then stopped while the

- 26 -

pipette resistance is measured and the motor position recorded ( $d=d_1$ , as shown in figures 1b and 3). If  $R_{tot}$  is outside a predetermined range, the experiment is aborted.

In step 04,  $R_{tot}$  is measured and if it is above a  
5 predetermined threshold, it is assumed that a cell is positioned on the pipette tip so suction is applied to the pipette and the logic controls the patch module motor to move the pipette tip further into the liquid within the capillary to a predetermined recording position ( $d=d_2$ , as  
10 shown in figures 1b and 3). The logic then moves to step 05 to test the G-seal.

If at the start of step 04,  $R_{tot}$  is less than the predetermined threshold, the logic assumes that there is no cell at the pipette tip. The logic then waits for a  
15 predetermined time interval  $t_1$  before controlling the patch module motor to move the capillary away from the pipette until  $R_{tot}$  is measured to be greater than the predetermined threshold, when the movement is stopped. It is believed that in this position the pipette tip is still  
20 in contact with the liquid in the capillary but only via a neck, or bridge, of liquid drawn out by surface tension between the capillary and the pipette. The logic then waits until  $R_{tot}$  drops to equal the predetermined threshold. The logic then controls the patch module motor  
25 to return the pipette tip to  $d=d_1$ , the position when it first contacted the capillary meniscus in step 03. If  $R_{tot}$  is then greater than the predetermined threshold it is assumed that contact with a cell has been made at the pipette tip, suction is applied to the pipette and the  
30 logic controls the patch module motor to move the pipette into the capillary to the predetermined recording position at  $d=d_2$ .

- 27 -

It is believed that waiting for the time interval  $t_1$ , which may be between 0.5 and 10 seconds, or preferably about 1 to 5 seconds, permits movement of the cells at the capillary tip, which is encouraged by the movement of the pipette tip to draw out the capillary meniscus.

If  $R_{tot}$  is still less than the predetermined threshold, the steps of waiting for time  $t_1$  and slightly moving the pipette are repeated for a predetermined number of iterations until a failure condition (step 20 ) is reached.

#### Control Logic - Third Embodiment

Figures 10 to 16 are flow charts illustrating a third embodiment of the control logic of the invention. Aspects of the third embodiment are, where appropriate, common to the first and second embodiments. The third embodiment, however, incorporates certain improvements resulting from experiments by the inventor.

Figure 10 is a flow chart showing all of the operations of the control logic, or software. This is termed the AutoPatch system. Figures 11 to 16, and figures 10a, 10b, 11a, 16a and 16b, are expanded flow charts for operations within the flow chart of Figure 10.

Figure 10 describes the setting up of the AutoPatch system, including the initialisation of all the relevant hardware. This involves the steps up to starting the test sweep 302, after which the steps of interface or cell location 304, Giga seal testing 306, whole cell detection 308, qualification 310 and experiment 312 are performed as described herein. Note that during the patching process, the movements of the capillary and the petri dish are locked together by the software in order to maintain their

- 28 -

positions relative to each other constant. The movement of the petri dish has no role in the patching process.

Initially the capillary and petri dish move towards the pipette at a rapid speed to a pre-set distant (step 314) in order to position the liquid/air interface within approximately 1mm of the pipette tip. This initial (rough) position is performed in order to shorten the time interval between starting the patching process and reaching the interface. The distance between the patch pipette tip and the liquid/air interface is initially larger than 1mm in order to allow sufficient space for loading of the pipette and the capillary. After rough positioning, the seal test pulse is started (step 316) and the translators (capillary and petri dish) are switched to slow speed prior to entry of the pipette into the liquid/air interface.

Figures 10a and 10b expand the safe quit routine 318 and the reset routine 320 used at various points by the control logic or software.

Figure 11 is a flow diagram expanding the interface or cell locator routine 304 of Figure 10. In this routine, suction is applied to the interior of the pipette, and the capillary and petri dish are moved towards the pipette at a slow rate (e.g. 10 micrometers/s) (step 322). The current is monitored after each sweep to determine when the interface is detected (step 324). The interface is detected by an offset in the baseline current or the appearance of a current pulse when the circuit is made 326. In contrast with conventional patch clamp, seal formation with the Interface Patch Clamp technique can occur virtually simultaneously with contact between the pipette tip and the external bathing solution at the

- 29 -

interface 328. The logic distinguishes between an open circuit (pipette tip in air) and a rapidly formed seal by monitoring the current trace for the appearance of a capacitance transient at the end of the seal test pulse (step 326). This transient is due to the pipette capacitance, which increases as the pipette is immersed into the solution at the interface. As the pipette tip touches or crosses the air/liquid interface the capacitance is likely to be undetectable amongst electrical noise. As the pipette tip moves into the liquid, however, the capacitance increases until it is detectable over the noise.

After the pipette has entered the interface, the capillary (and movement of the petri dish) is stopped, a junction null may be performed and the pipette resistance is monitored. The presence of a cell on the pipette tip is indicated by an increase in pipette resistance which must reach a pre-set value before the logic will allow the system to proceed to the next step. After a cell has been detected the Giga seal test 306 is undertaken.

Figure 11a expands the interface hit routine 330 of the interface or cell location routine 304 of Figure 11.

Figure 12 is a flow diagram of the Giga seal test routine 305 carried out after contact with a cell has been detected in the interface or cell location routine 304. The Giga seal test routine 306 comprises a repetitive loop 332 in which the level of suction applied to the pipette interior is increased in pre-set increments 334 and times 336 while monitoring the change in pipette resistance. The suction is increased until the maximum vacuum is obtained or Giga seal formation occurs. The loop ends and the suction is switched off 338 if either of these

- 30 -

conditions is satisfied. If the maximum suction has been applied but a Giga seal has not formed the loop is repeated until Giga seal formation or time-out. The formation of a Giga ohm seal is required in order to allow  
5 progression to the next stage.

Figure 13 expands the whole cell detection routine 308 of Figure 10. In this routine, the transients due to the pipette capacitance are cancelled and this is followed by a repetitive loop 340 in which the level of suction  
10 applied to the pipette interior is increased in pre-set increments 342 and times 344 while monitoring the current 346 for the appearance of capacitance transients. These transients are due to the charge and discharge of the cell membrane capacitance and are indicative of the formation  
15 of whole cell recording mode.

Figure 14 expands the qualification routine 310 of Figure 10. In order to qualify for use in an experiment the cell must exhibit a voltage gated (or other) current equal to or greater than a pre-determined amplitude and polarity  
20 348 in response to a test pulse. Qualification proceeds until the cell qualifies or timeout is hit. During qualification the quality monitor 350 is run also.

Figure 15 is a flow diagram of the quality monitor routine 350 of Figure 14 and 16, described below. It comprises a  
25 repetitive loop in which the pipette suction is varied in response to measurements of series resistance (RSeries) and current (IMon). Current flow through the cell membrane via the pipette generates a voltage error due to RSeries. The value of RSeries often increases during  
30 whole cell recording and this effect can be reduced by the application of suction to the pipette interior 352. An increase in the value of the current at the holding



- 31 -

potential (usually - 60mV) indicates loss of the Giga seal and this can be caused by excessive suction. Acceptable values for RSeries and IMon are entered in the settings for the software. The quality monitor runs both during  
5 the qualification stage and the experiment.

Figure 16 expands the experiment routine 312 of Figure 10. In the experiment routine, movements of the petri dish and capillary required to carry out drug application by the method shown in figures 4-7 inclusive are carried out.  
10 During these movements the dish is filled with external solution/external solution plus drug via solenoid operated flow control values 358. Before drug is added, the current evoked by the test pulse (or pulses) must be reproducible within a pre-determined percentage (entered  
15 in the settings).

Figures 16a and 16b expand the trace run routine 354 and the stabilise routine 356 of the experiment routine 312.

- 32 -

Figure 17 shows recordings of current (1) and voltage (2) obtained from an automated patch clamp system (AutoPatch) showing the formation of a Giga Seal under software control using the Interface Patch Clamp technique.

5 Recordings from an MK1 cell.

Figure 18 shows recordings of current (1) and voltage (2) from an automated interface patch clamp system (AutoPatch) showing the increase in capacitance transient observed after moving the capillary to position d2 (records a and  
10 b). (c) and (d) were obtained after automatic compensation for pipette capacitance and a change in the holding potential to -60mV. Recording obtained from the same cell as Figure 17.

Figure 19 shows recordings of current (1) and voltage (2) from an automated interface patch clamp system  
15 (AutoPatch). Records (a) and (b) were obtained after Giga seal formation (cell attached patch mode) using the interface patch clamp technique. The application of suction to the pipette interior under full software  
20 control ruptured the membrane patch to obtain the whole cell recordings shown in (c) and (d). The establishment of the whole cell mode of recording is shown by the presence of the large capacitance transients on the current trace. Recording obtained from the same cell as  
25 figures 17 and 18.

Figures 20 shows recordings of membrane current (1) and voltage (2) in whole cell recording mode obtained using an automated patch clamp system (AutoPatch) which employs the interface patch clamp technique. The holding potential  
30 was -60mV and (a) and (b) show outward potassium currents (Kv1.1) in response to a voltage step to +30mV. Recording obtained from the same cell as figures 17, 18 and 19.

- 33 -

Figure 21 shows the effect of the potassium channel blocking drug tetraethylammonium (TEA) on the potassium current recorded from MK1 cell in whole cell recording mode obtained using the Interface Patch Clamp technique. After the establishment of whole cell recording mode the cell was positioned in a recording dish by the method shown in figures 4 to 7 and described in the text.

- (a) shows the current obtained in normal extracellular solution
- (b) shows the effect of replacing the solution in the dish with extracellular solution containing TEA (5mM).
- (c) the blocking effect of TEA was reversed by washing.

It will be readily appreciated by those skilled in the art that:

1. The stability of recording using the interface patch clamp technique may be superior to that of conventional patch clamping. The greater stability of interface patch clamping is because the cell is held by the patch pipette alone. In conventional patch clamp recordings the cell is held by the path pipette and a solid substrate and vibration tends to move the pipette relative to the substrate causing loss of the G-seal. The interface patch clamp is, in contrast to conventional patch clamp apparatus, relatively insensitive to vibration during drug application.
2. This method of drug application could be applied to a plurality of recording pipettes/capillaries and form the basis for a high throughput electrophysiological assay system. It will readily be appreciated that the Interface Patch Clamp technique could be used with multiple pipettes and multiple capillaries in a

- 34 -

manner in which each pipette enters its respective aligned capillary either individually in sequence or all together. Although not currently preferred, a single pipette could be used which is caused to enter more than one capillary sequentially. Multiple patch clamp recordings could be made either sequentially or simultaneously, depending on the application.

As was mentioned above, it is not essential to the general principle of the invention to use a capillary in order to create a column of liquid which gives rise to a liquid/air interface at which cells can be located. Other ways can be envisaged in which the same effect can be achieved. For example, as shown in figure 22, a droplet or "blob" of liquid may be provided on a support surface. The surface has a hole through it and the droplet covers the hole. Surface tension prevents the liquid from the droplet dropping through the hole. Within the droplet cells are suspended. This allows access to the droplet and the cells contained therein by a suitable electrode such as a patch pipette. In the arrangement shown in figure 22, means are provided for flow of other liquids in to and out of a dish or other container of which the support surface with the hole in it forms a wall. Once a cell has been attached to the electrode, other liquids may be introduced into the container either in batch mode or in flow-through mode in order to result in the cell being exposed at its external surface to the surrounding liquid. Clearly in such an arrangement, the original liquid and the remaining un-attached cells will tend to be washed away from the area of the electrode/pipette.

- 35 -

It is within the scope of the invention that droplets might be provided on non-perforate support surfaces. As shown in figure 23, the effect of surface tension may be to allow droplets of a suitable liquid to adhere  
5 automatically to the underside of a suitable support surface. The support surface might for example be a cover slip of glass or other material. Droplets in which cells are suspended provide the air/liquid interface according to the invention and consequently may be used in a method  
10 of interface patch clamping as described above.

As has already been mentioned, the arrangement shown in Figure 23 as well as the arrangement shown in Figure 24 allows for the formation of a matrix of cell suspensions so that multiple electrodes can be multiplexed to take  
15 readings either simultaneously or sequentially (as well as singly).

It will be appreciated by those skilled in the art that a conventional glass "patch pipette" could be replaced by an equivalent electrode. It is considered to be within the  
20 scope of the present invention that the electrode might be either a single region or a matrix of regions on a sheet of material (such as a silicon wafer) which incorporates a microstructure to which a cell can be attached and which would provide the necessary electrical connection. For  
25 example, as shown in Figure 25, microstructures could be etched on to a silicon wafer (e.g. an oxidised silicon wafer), which microstructures would be designed and adapted to be able to capture a cell from the liquid/air interface of an arrangement according to the present  
30 invention. Thus, the performance and advantage of the invention is not limited to the currently preferred

- 36 -

conventional glass patch pipette but would include functionally equivalent means.

As has been described before, a drug in liquid solution can be applied to the cell in a number of ways. For example the drug could be applied via the capillary if the air interface is formed in a capillary tube. Alternatively the drug can be applied by perfusion into a dish (as described with reference to figures 4 to 7). Furthermore, perfusion could be achieved by flowing the drug-containing liquid through a dish or container as shown in Figure 22.

A further arrangement for drug application is shown in Figure 26a and 26b. In this case the electrode (for example the patch pipette) penetrates through the lower wall of a well. A suspension of cells is loaded in to a capillary tube as previously described. Attachment of a single cell to each pipette tip follows, as described before. Once cells are attached to the pipette tips the capillary tubes containing the remainder of the cells in suspension can be removed. Subsequently, a drug solution (23) is dispensed into each well (Figure 26b) and patch clamp measurements can then be carried out on the cell in the environment of the surrounding drug solution.

#### Optimisation of Patch Clamping Conditions

Those skilled in the art will appreciate that within the general teaching contained herein for the interface patch clamping method and apparatus, it may be necessary to optimise certain conditions for patch clamp measurements. For example the concentration and packing density of cells in the suspension may need to be optimised. Furthermore, the cells and/or solutions may be temperature sensitive

- 37 -

and an optimum temperature of operation may need to be determined. Since the invention relies on the formation of a liquid/air interface at which the cells are located, it may be necessary to optimise the osmolarity of the  
5 suspending liquid medium in order to achieve the optimum level of surface tension etc.

- 38 -

**Claims**

1. A method for providing a cell attached to a patch clamp electrode and having a high resistance (Giga Ohm) electrical seal between an area of the cell membrane and the electrode, which includes the steps of:
- 5 i) providing a suspension of cells in a liquid;
  - ii) causing the formation of a layer of cells at the interface between the air and the liquid in which the cells are suspended;
  - 10 iii) bringing the patch clamp electrode into contact with the interface by moving one or both of the electrode and the interface respectively together;
  - iv) contacting the electrode with a cell in the cell layer at or near the interface; and
  - 15 v) causing attachment of the cell to the electrode.
2. A method according to claim 1 for providing a cell attached to the tip of a patch clamp pipette and having a high resistance (Giga Ohm) electrical seal between an area of the cell membrane and the tip, which includes the steps of:
- 20 i) providing a capillary tube containing a suspension of cells in a liquid;
  - ii) causing the formation of a layer of cells at one end of the capillary tube at the interface between the air and the liquid in which the cells are suspended;
  - 25 iii) bringing the tip of the patch clamp pipette into contact with the interface by moving one or both of the pipette and the tube respectively together along a common axis of movement;
  - 30 iv) contacting the tip with a cell in the cell layer at or near the interface; and
  - v) causing attachment of the cell to the tip.



- 39 -

3. A method according to claim 1 or 2 in which the liquid in which the cells are suspended is an extracellular physiological solution.
4. A method according to claim 1 or 2 in which the layer  
5 of cells is several cells deep and loosely packed.
5. A method according to claim 2 in which the layer of cells is formed by mounting the capillary tube in an essentially upright orientation and allowing the suspended cells to sediment to the downward end of the tube to  
10 collect there in a layer.
6. A method according to claim 2 in which the capillary tube is mounted essentially upright with the interface at a lower open end of the tube and the pipette is mounted essentially upright with the tip upwardly pointing.
- 15 7. A method according to claim 2 in which the capillary tube and pipette are concentrically mounted with the capillary tube in a fixed position and the pipette movable along the common axis.
8. A method according to claim 2 in which the capillary  
20 tube and pipette are concentrically mounted with the pipette in a fixed position and the capillary tube movable along the common axis.
9. A method according to claim 2 wherein gentle suction is applied to the pipette during contact with the  
25 interface and during the step of contacting the tip with a cell.
10. A method according to any preceding claim, in which

- 40 -

contact between the pipette tip and the air/liquid interface and/or subsequent movement of the pipette tip into the liquid is detected by monitoring pipette capacitance.

5      11. A method according to any preceding claim, in which  
if no cell is contacted at or near the interface as or  
within a predetermined time after contact between the  
pipette and the interface, the pipette is withdrawn from  
the interface and then moved back to the interface to  
10      repeat the attempt to contact a cell.

12. An apparatus for carrying out the method of any  
preceding claim which is a computer controlled apparatus  
including the following elements:  
i) a patch clamp amplifier;  
15      ii) a source of variable suction for a patch clamp pipette  
under the control of the patch clamp amplifier;  
iii) a holder for a capillary tube to be mounted  
vertically;  
iv) a holder for a patch clamp pipette to be mounted  
20      vertically in the same axis as the capillary tube in an  
inverted orientation with the tip pointing upwardly;  
v) a manipulator for controlling relative movement of the  
capillary tube and pipette along a common axis of movement  
under feedback control from the patch clamp amplifier and  
25      allowing for the tip of the pipette to enter a downwardly  
facing end of the capillary tube.

13. An apparatus according to claim 12 which includes an  
array of a multiplicity of capillary tubes and an array of  
a multiplicity of pipettes.

- 41 -

14. An apparatus according to claim 12 or 13, comprising a pipette capacitance sensor for sensing pipette capacitance as the pipette tip contacts an air/liquid interface at the end of the capillary tube and enters the liquid in the capillary tube during operation of the apparatus.
15. A computer-program-controlled patch clamping process for carrying out the method of any of claims 1 to 11.
16. A computer-program-controlled patch clamping process for controlling the apparatus of claim 12, 13 or 14.
17. A computer-readable medium carrying a computer program for controlling a computer to implement the method any of claims 1 to 11 or to control the apparatus of claim 12, 13 or 14.
18. A method for controlling a computer by means of a computer program for implementing the method of any of claims 1 to 11.

1 / 33

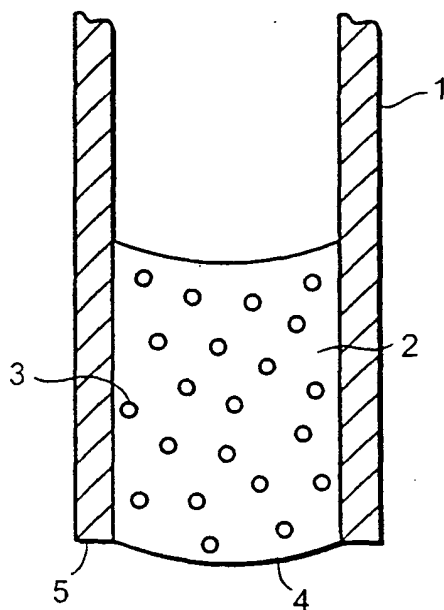


FIG. 1a

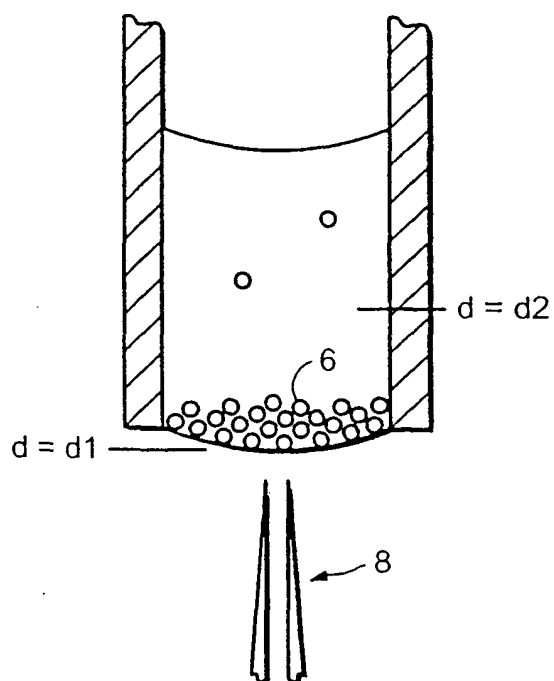


FIG. 1b

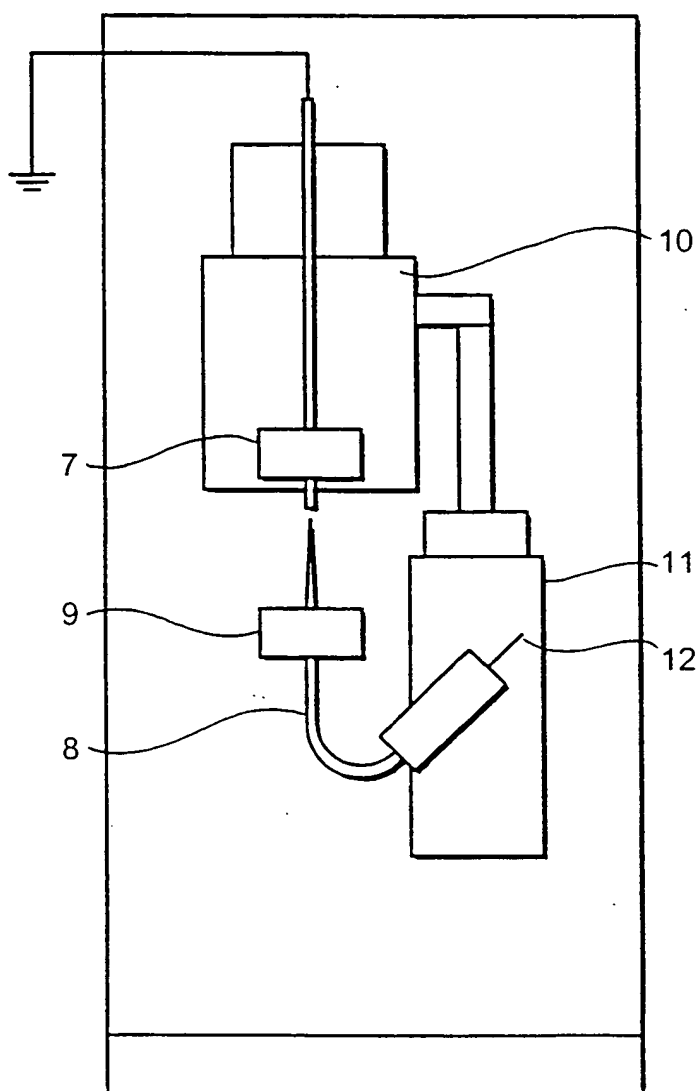


FIG. 2

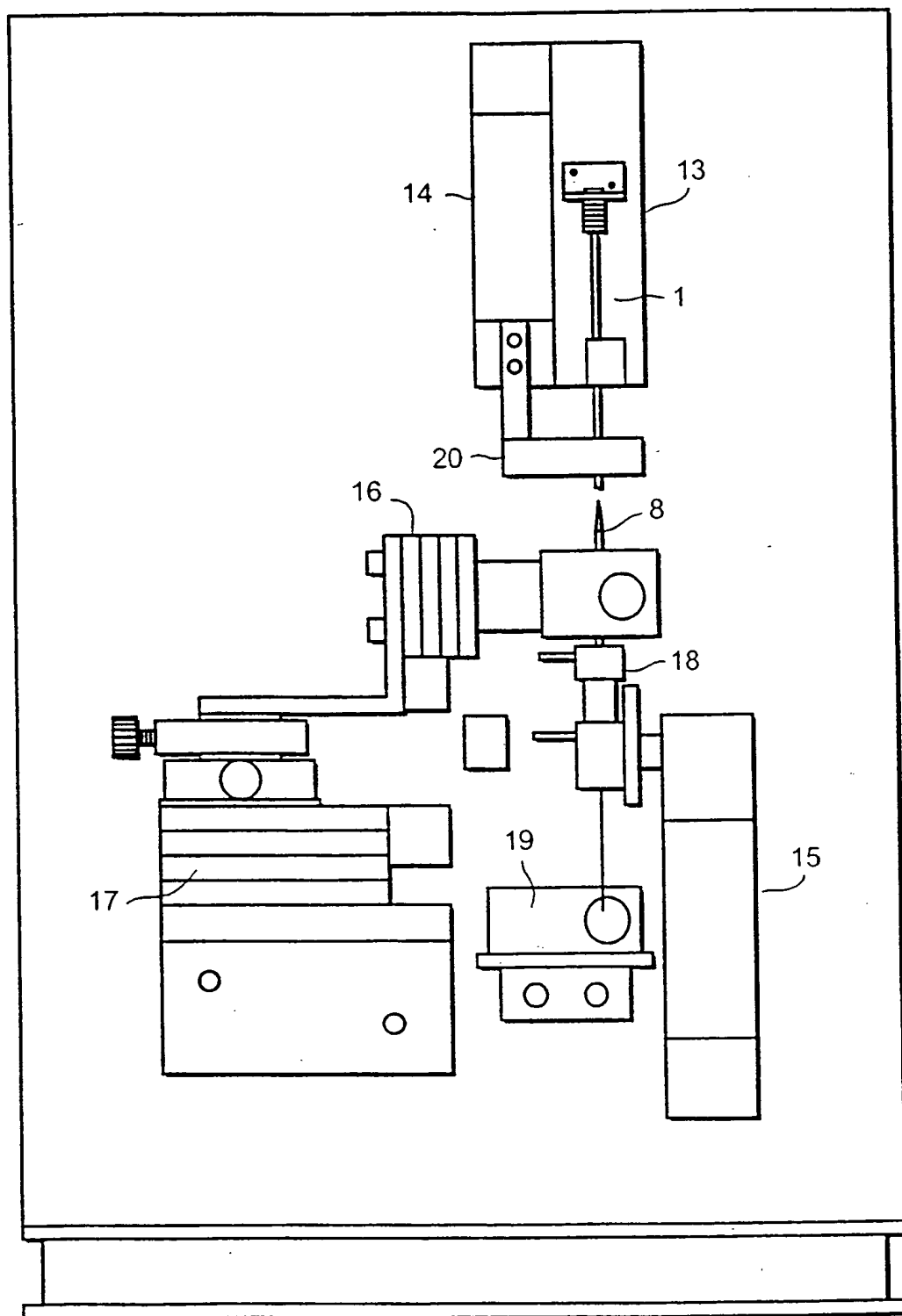


FIG. 2a

4 / 33

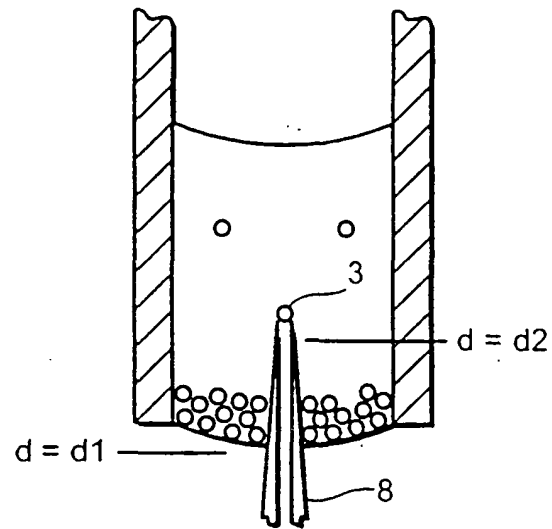


FIG. 3

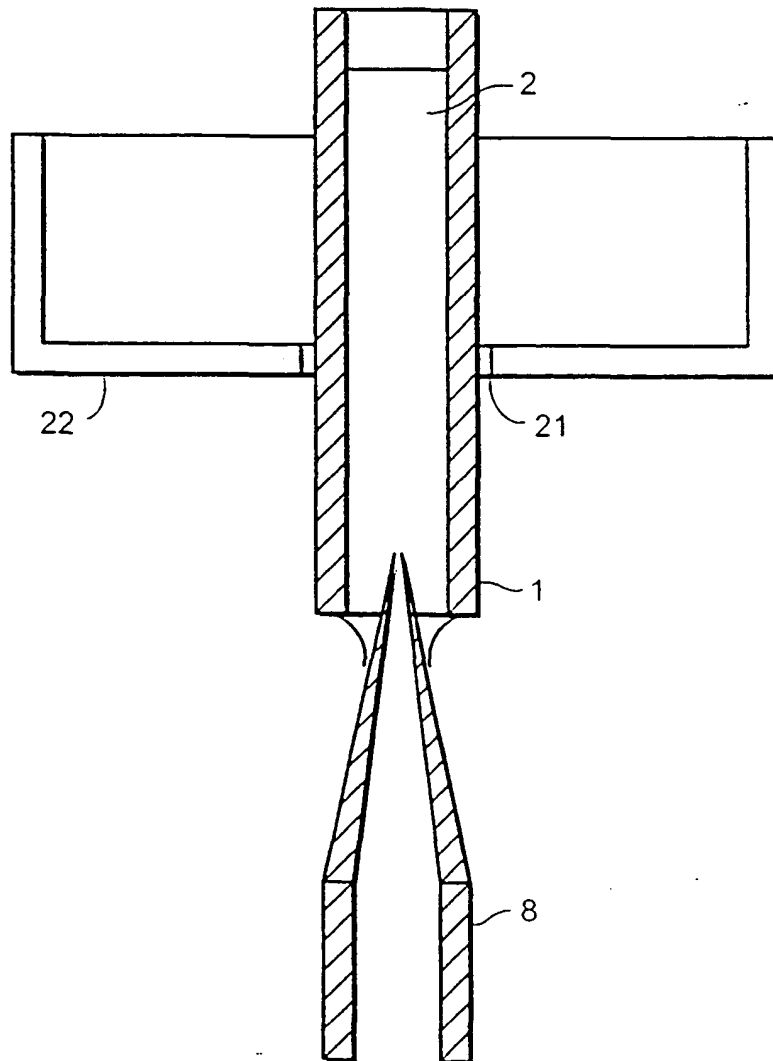


FIG. 4

5 / 33

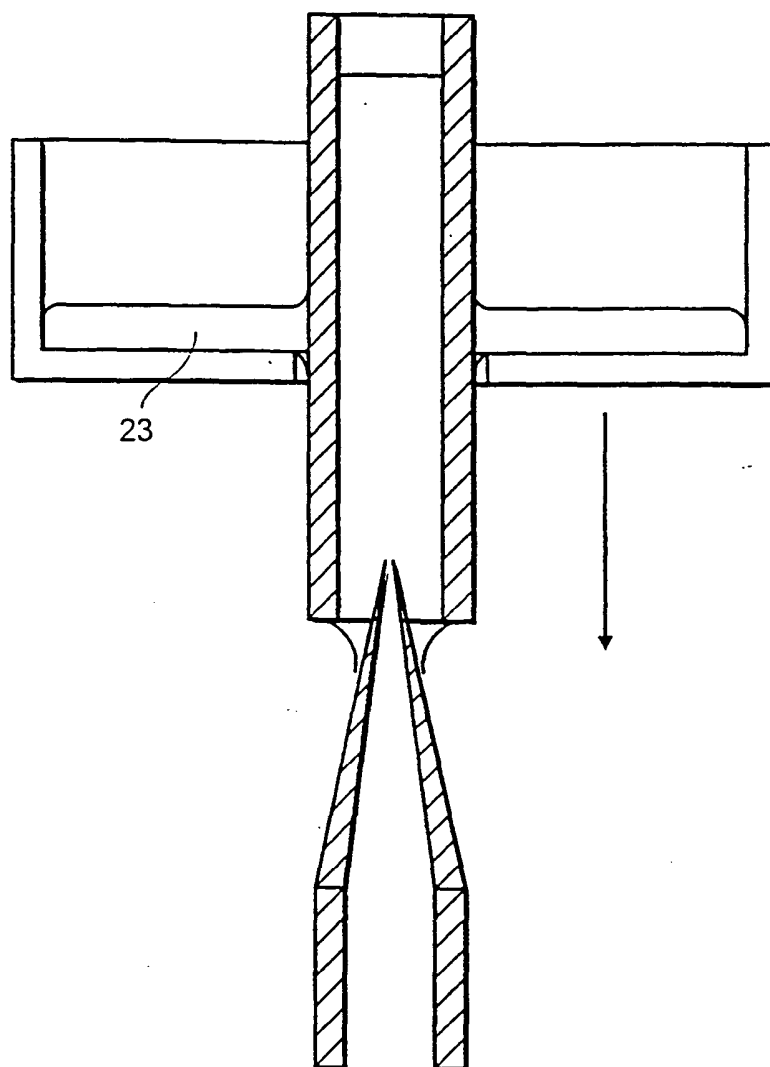


FIG. 5



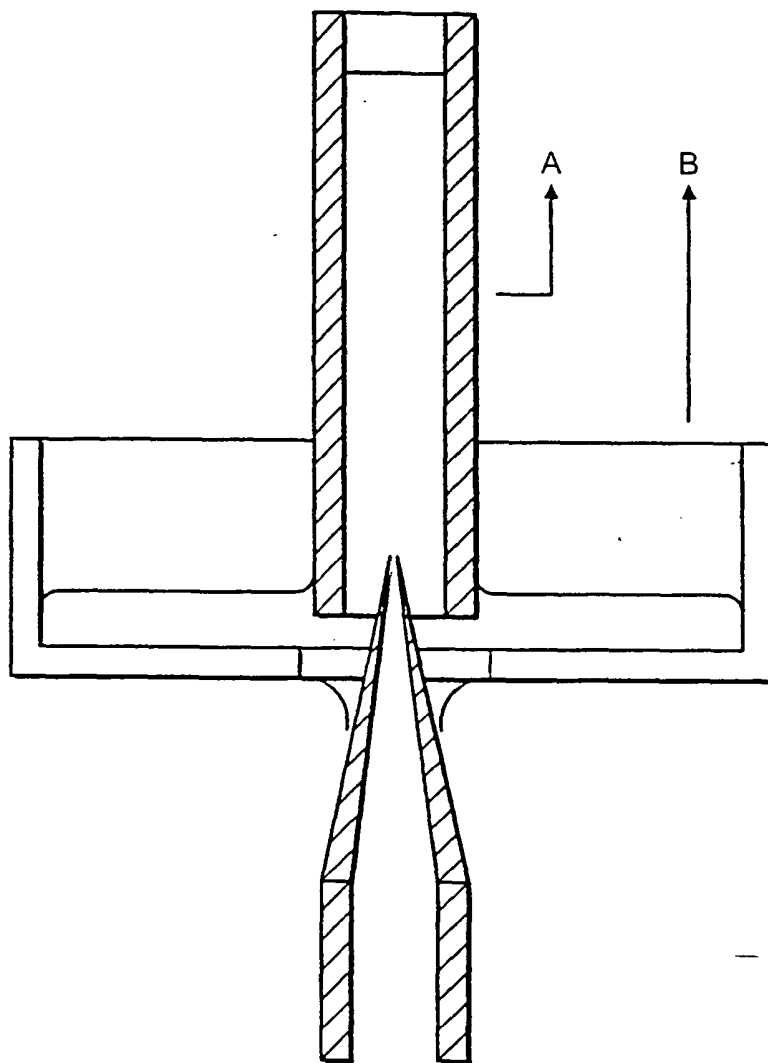


FIG. 6



8 / 33

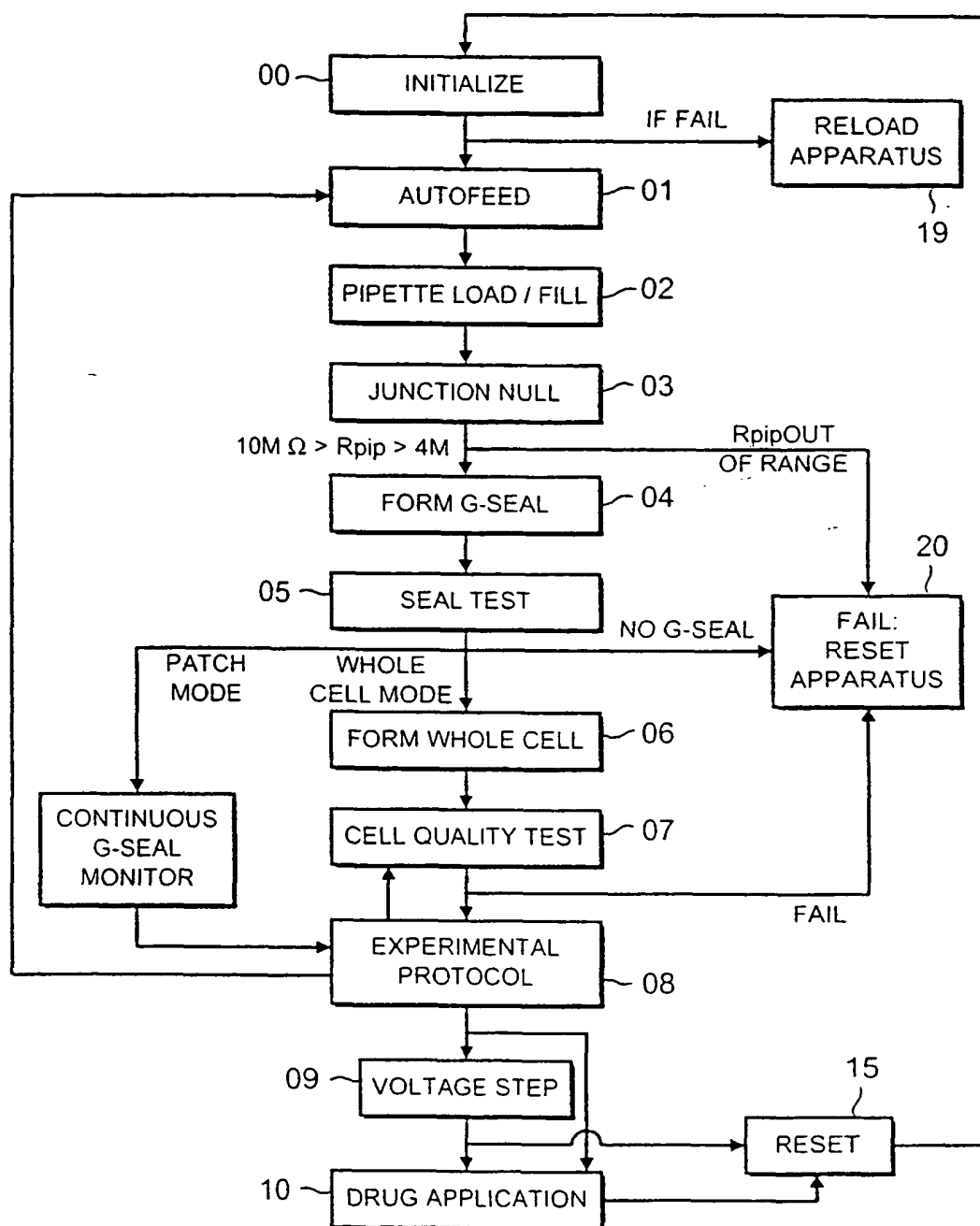


FIG. 8

9 / 33

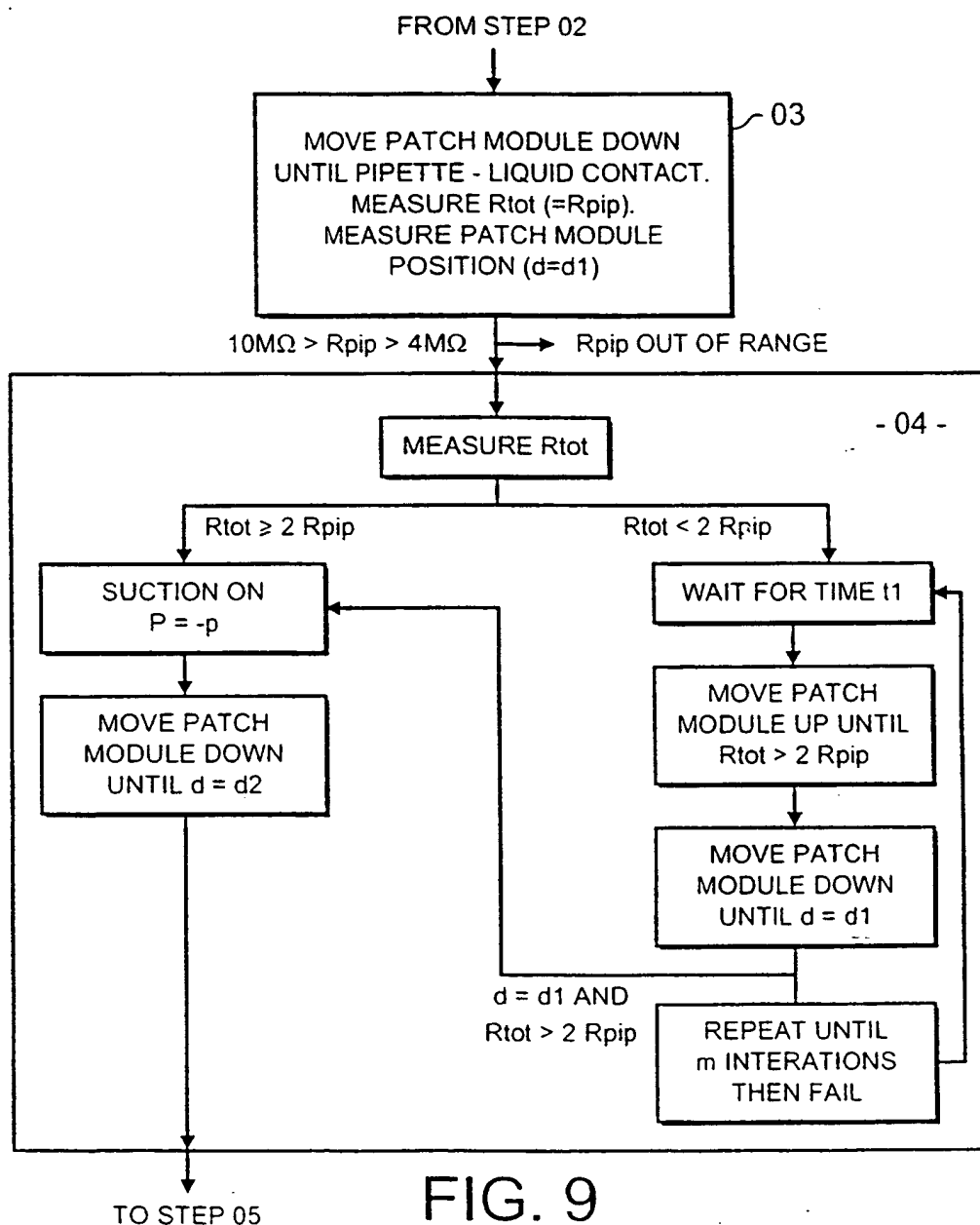


FIG. 9

FIG. 10

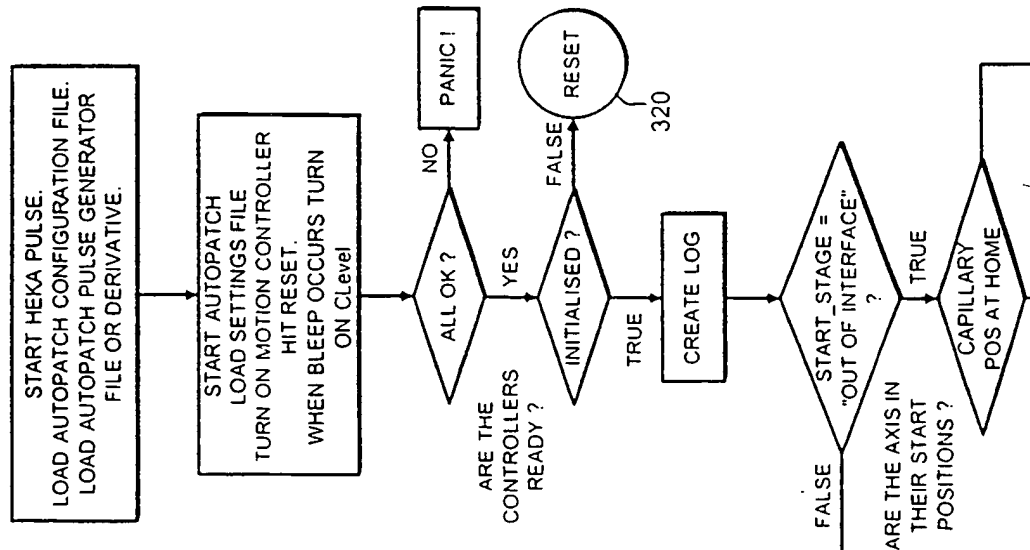


FIG. 10a

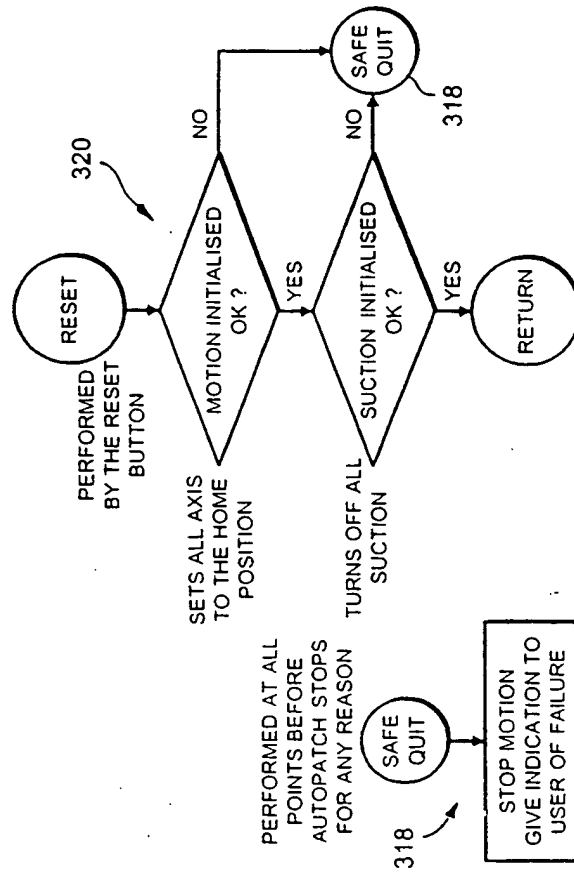


FIG. 10b

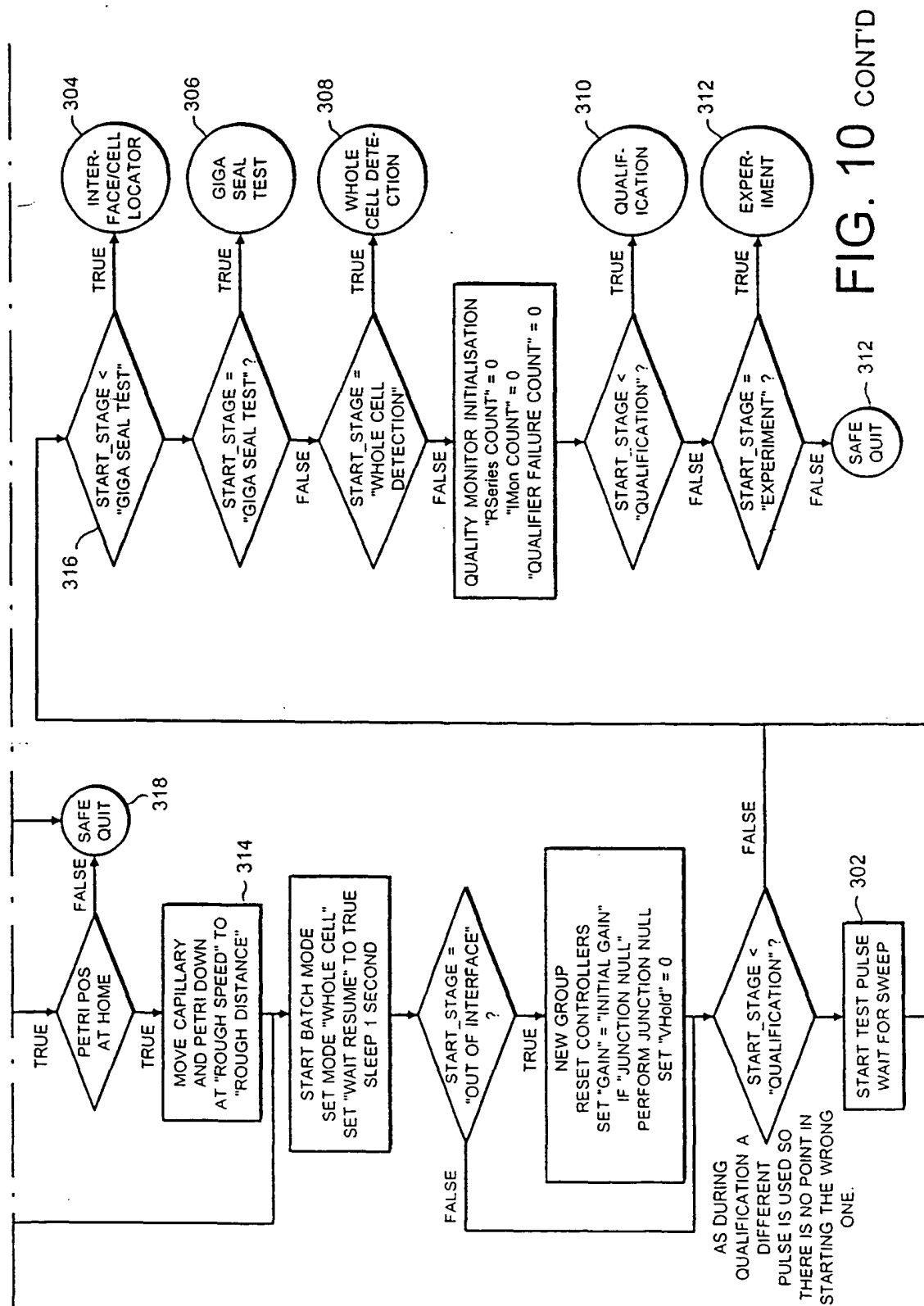
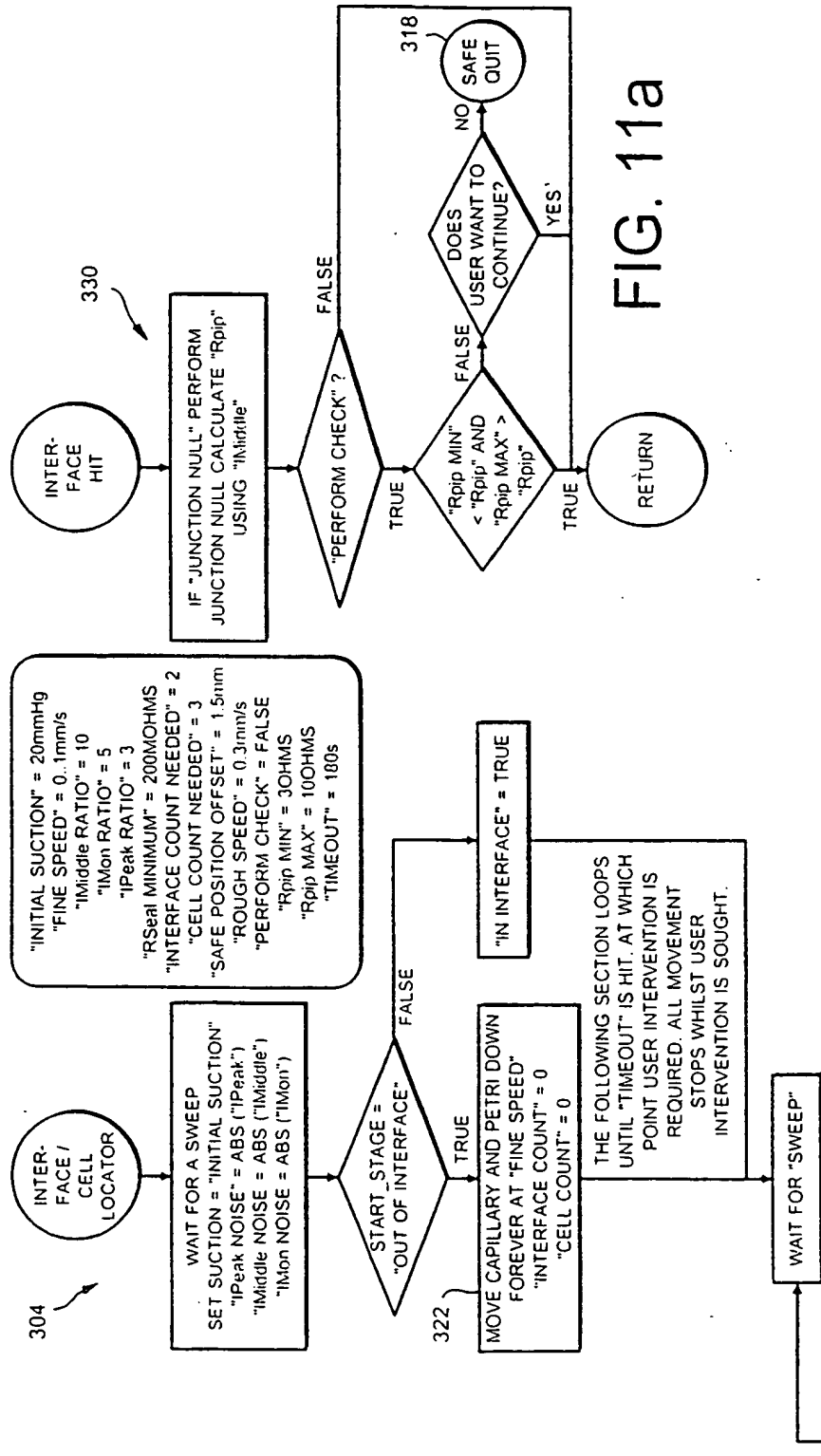


FIG. 10 CONT'D

FIG. 11





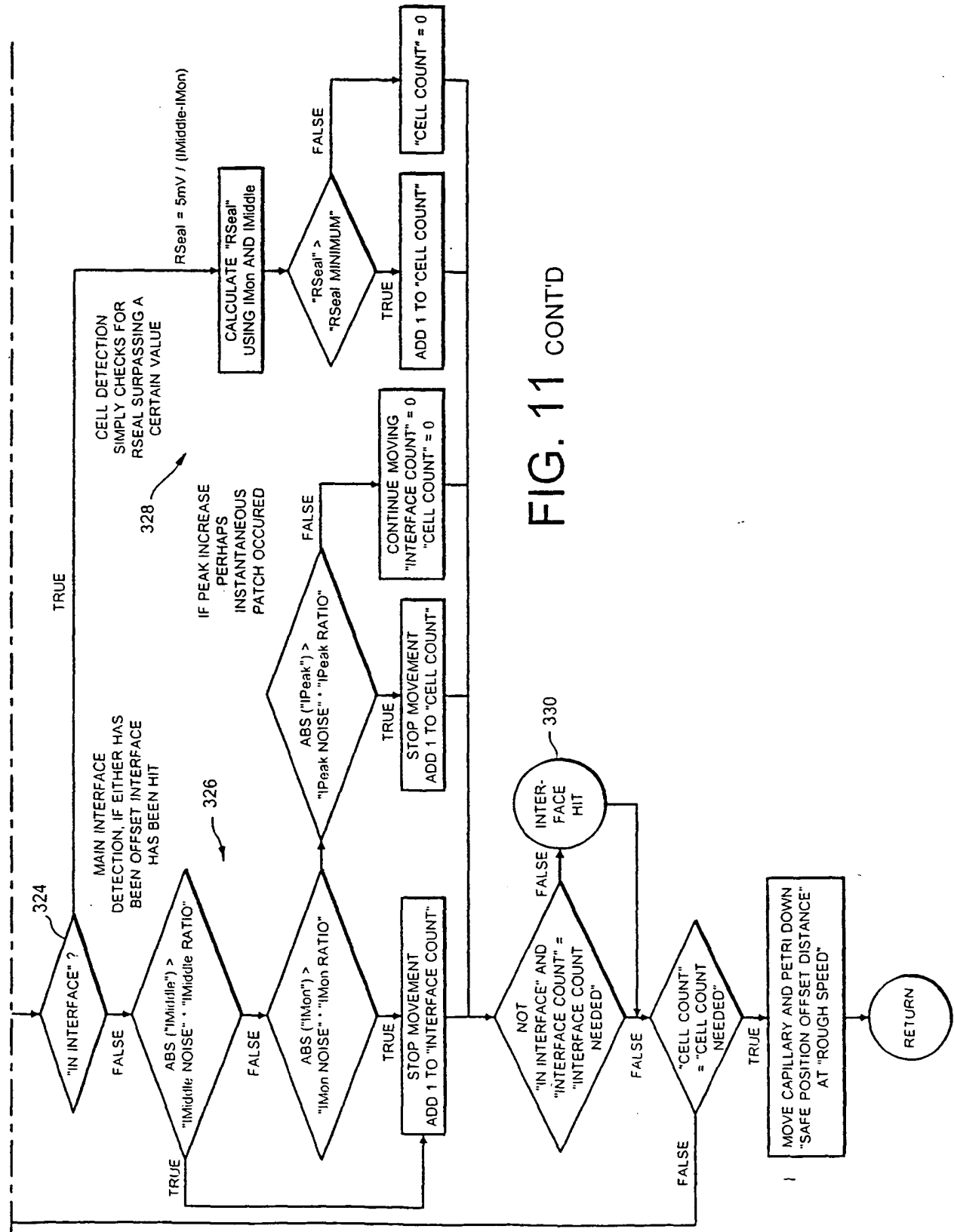



FIG. 11 CONT'D




FIG. 12

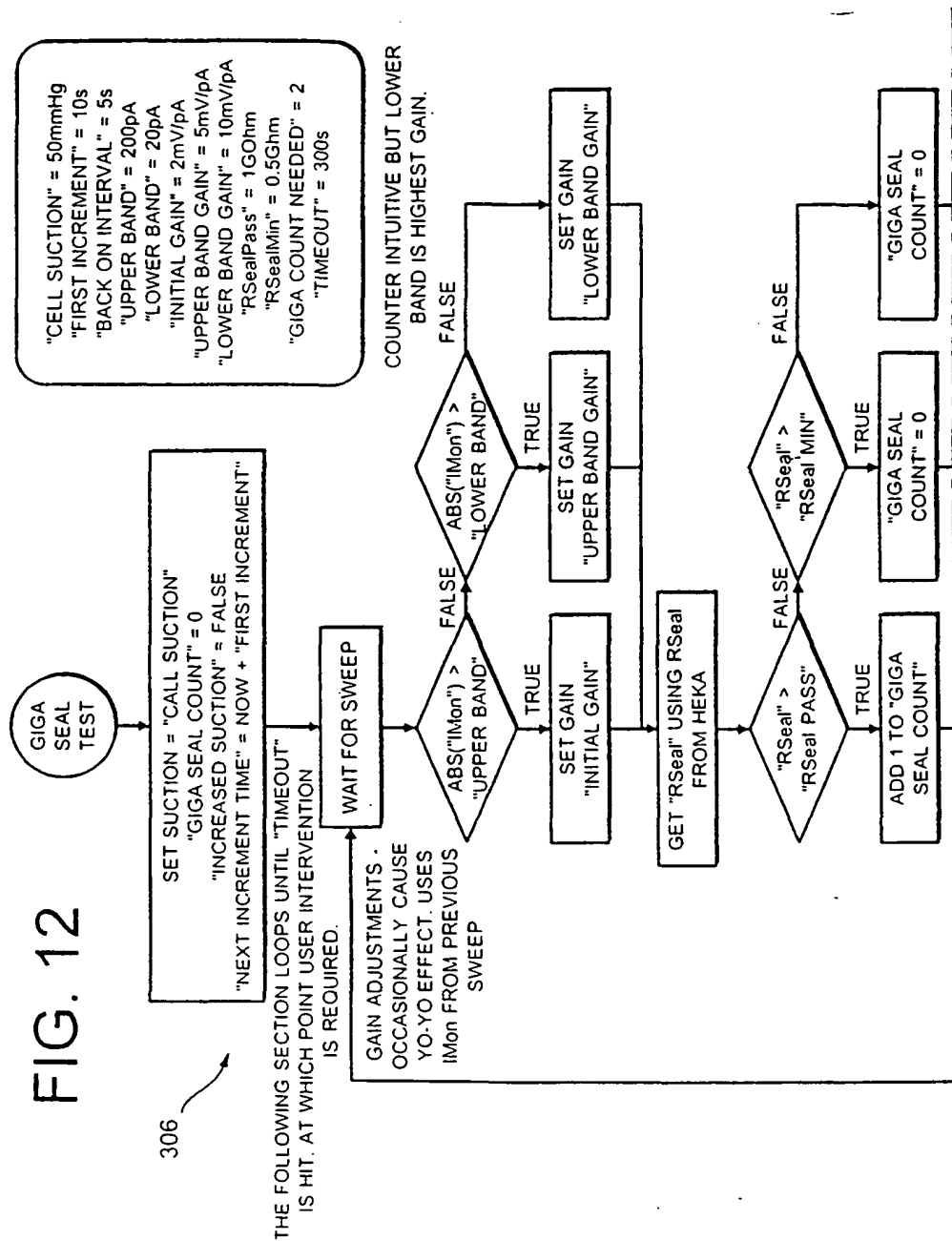


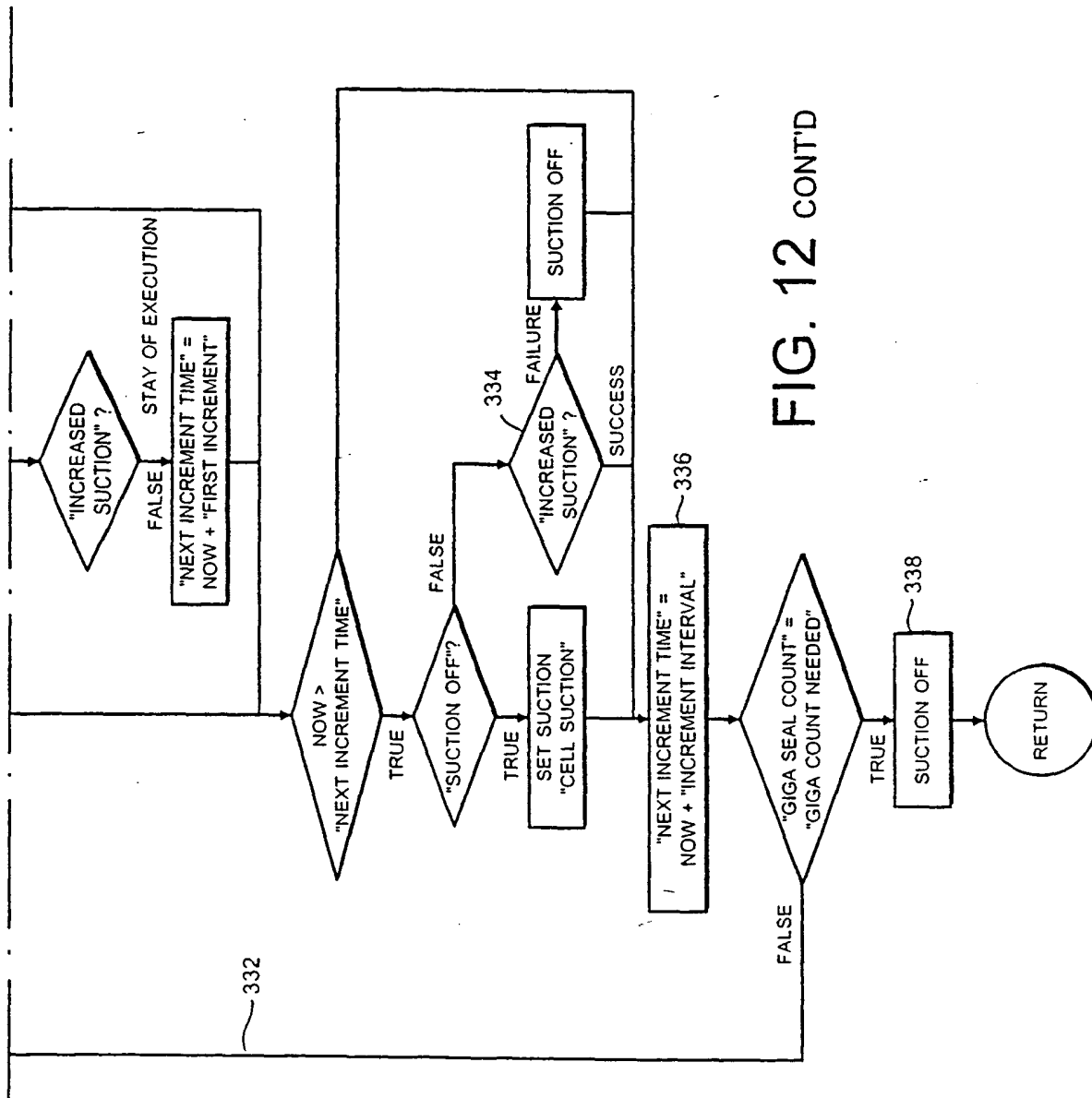


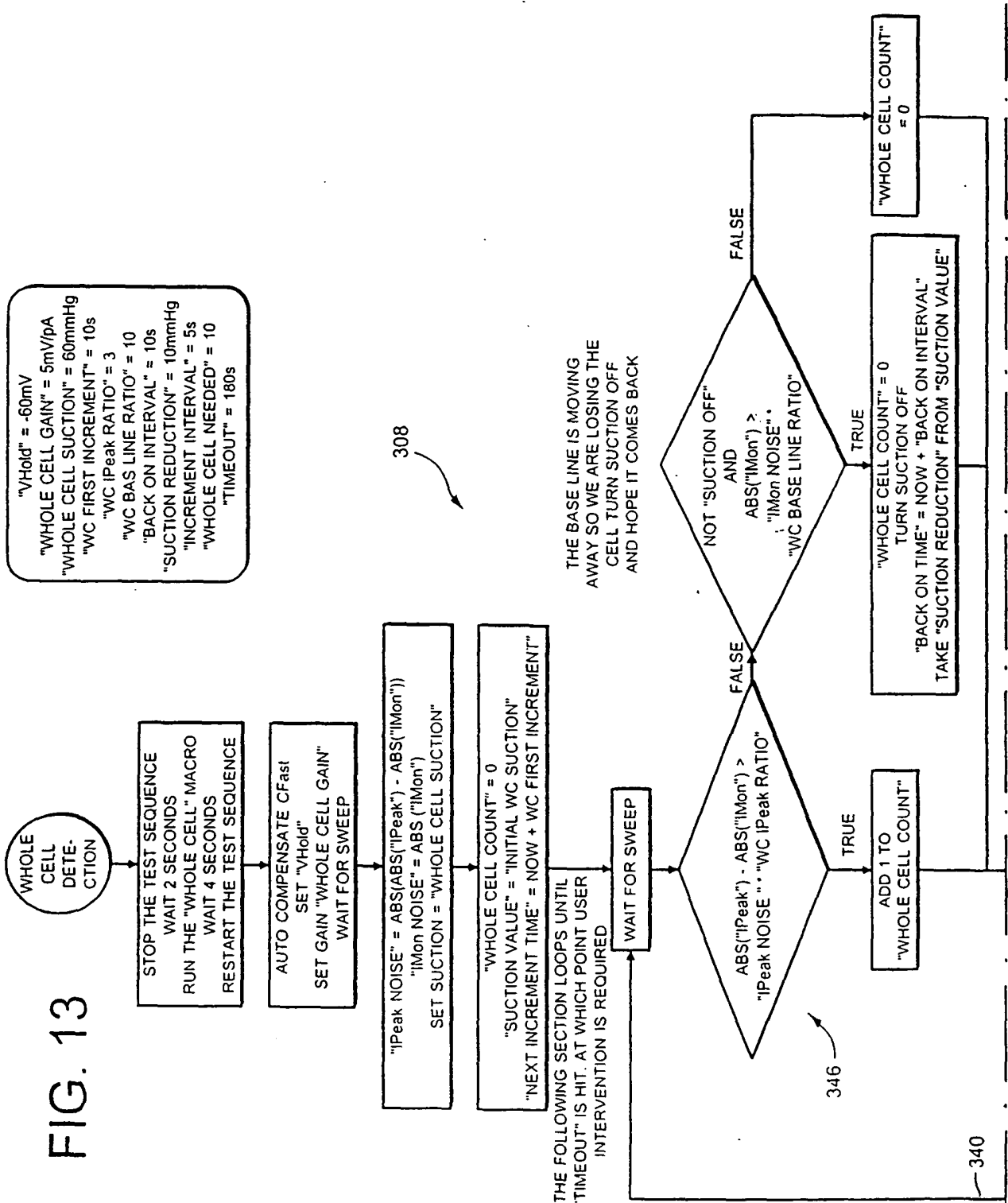


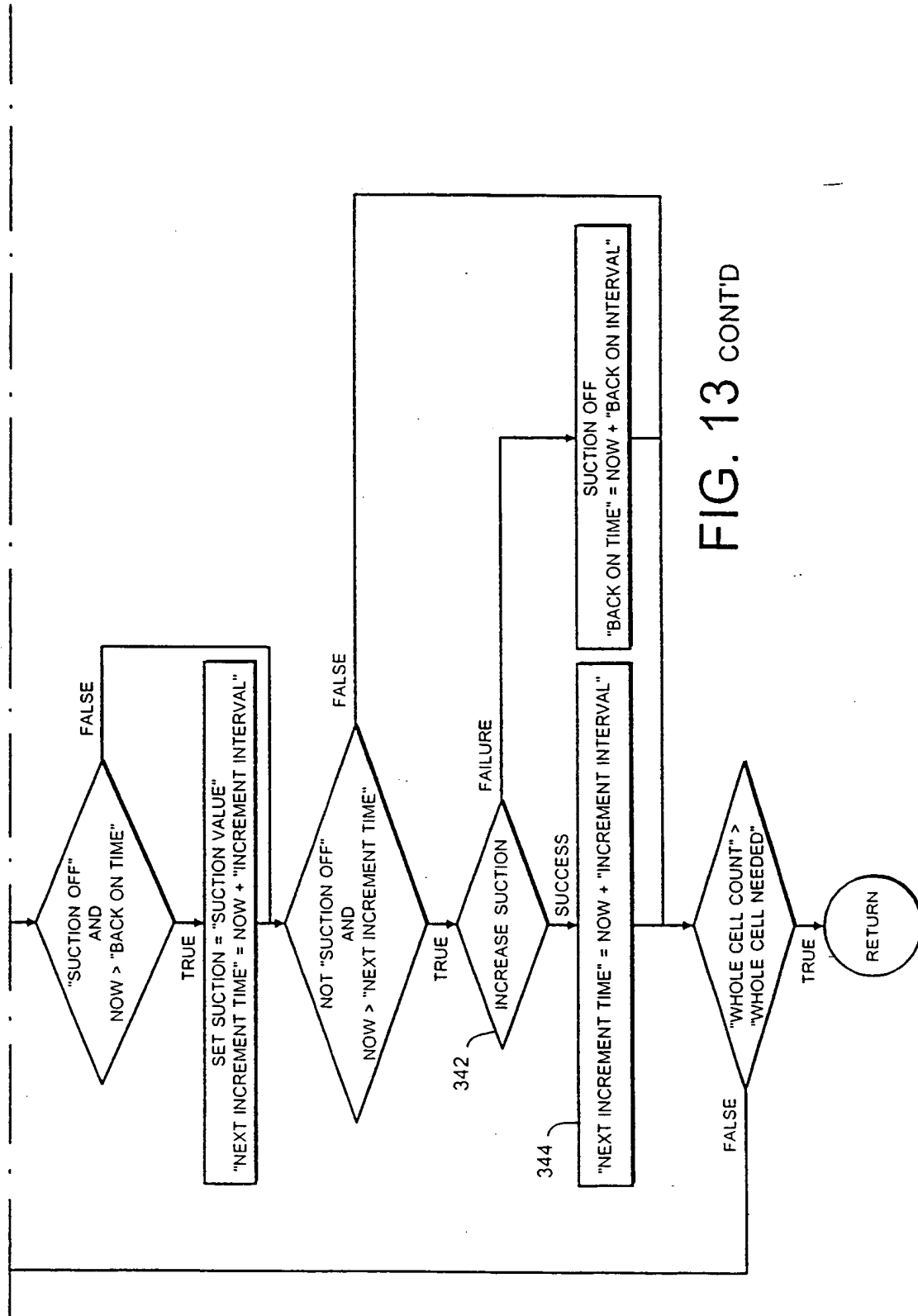

18 / 33

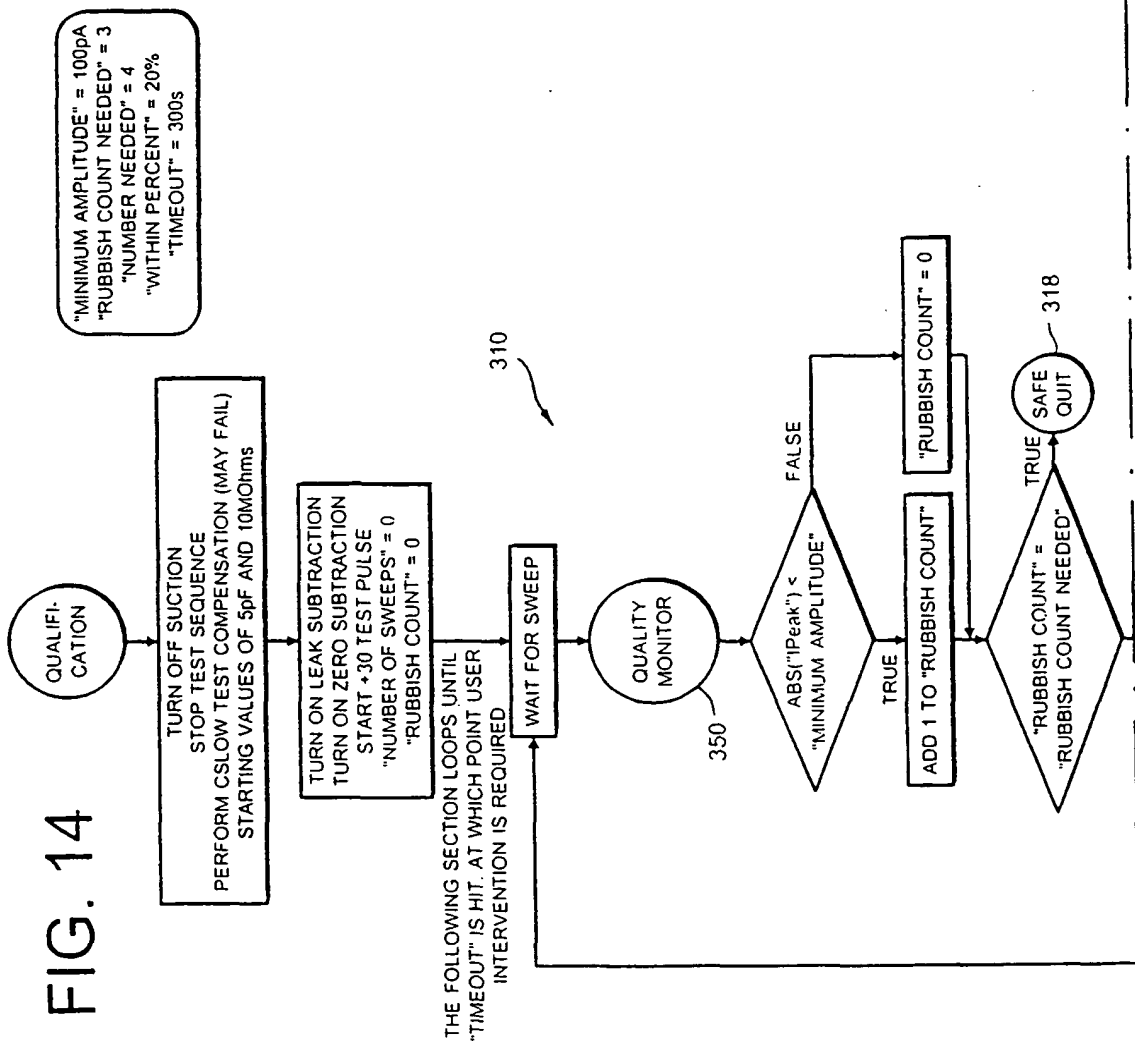

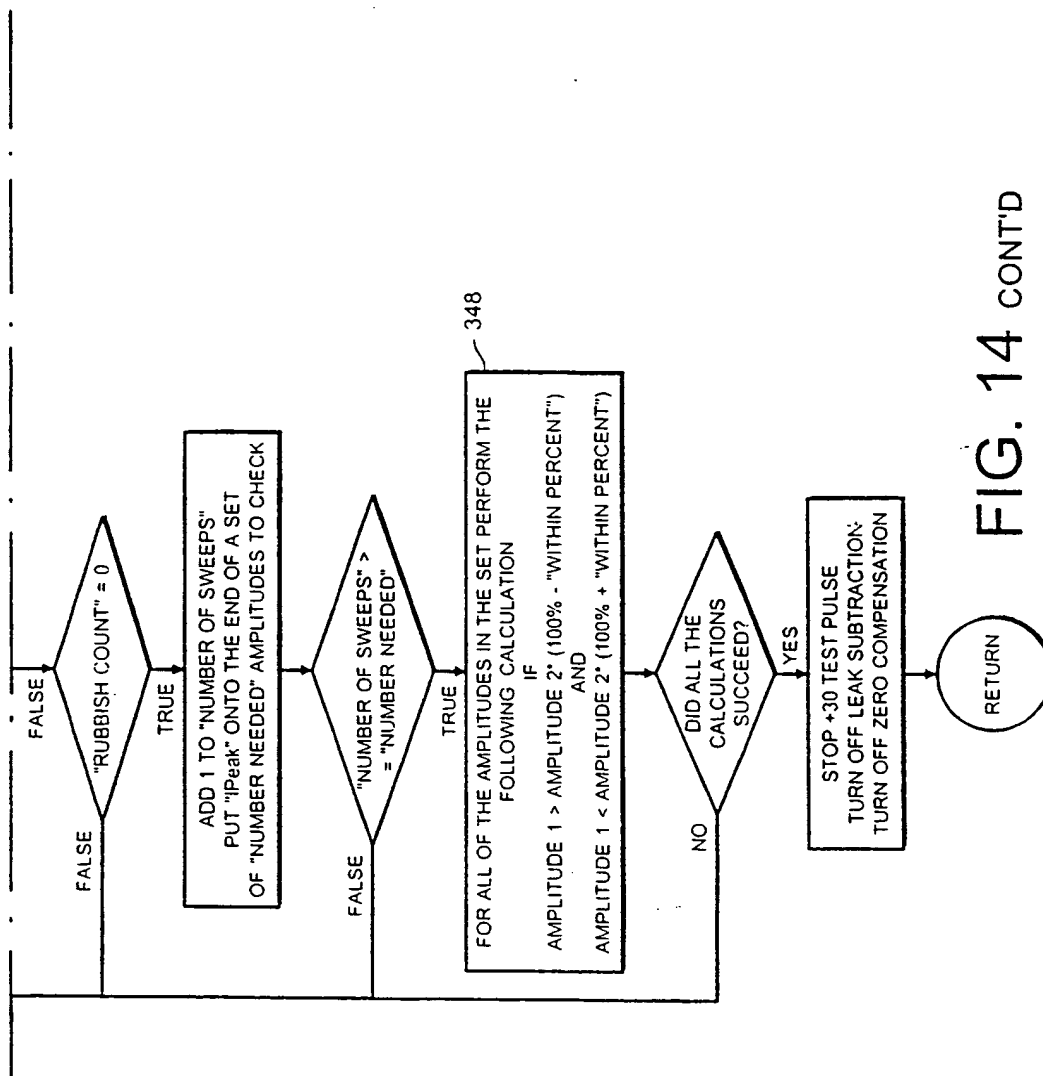

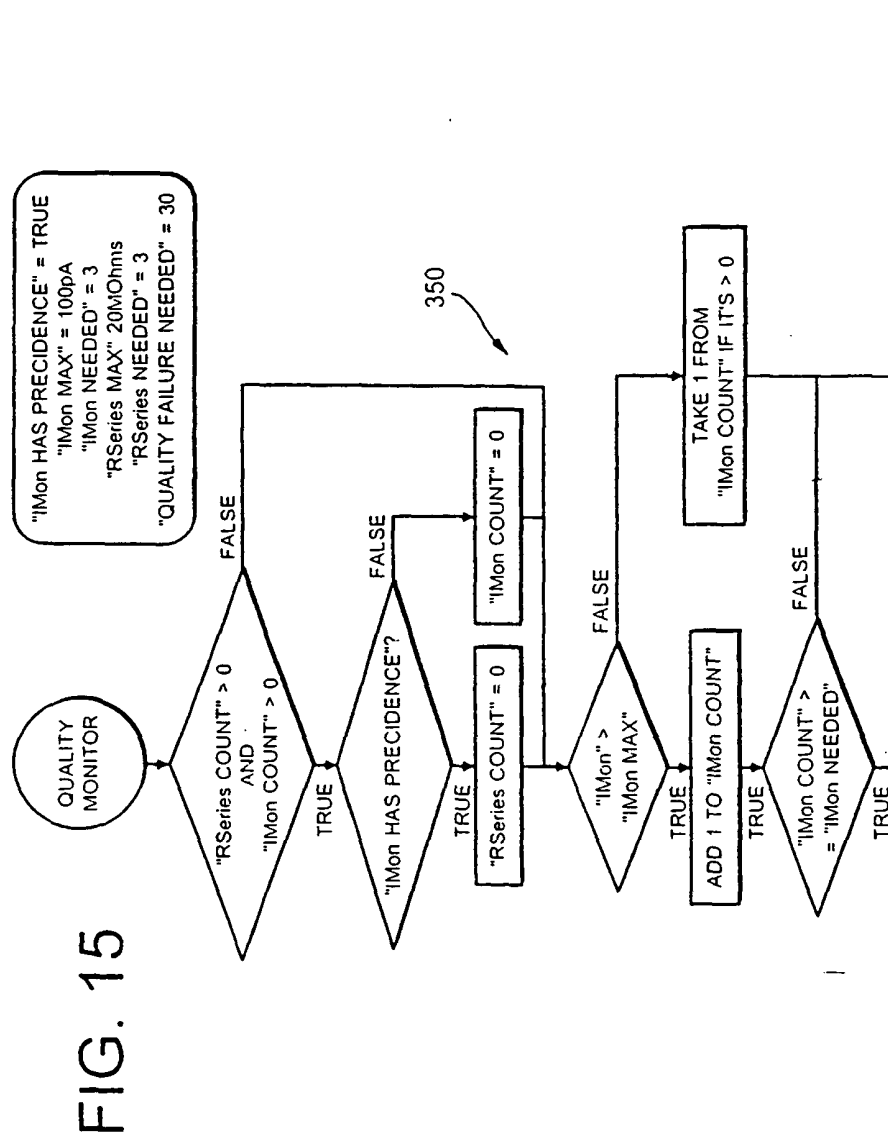



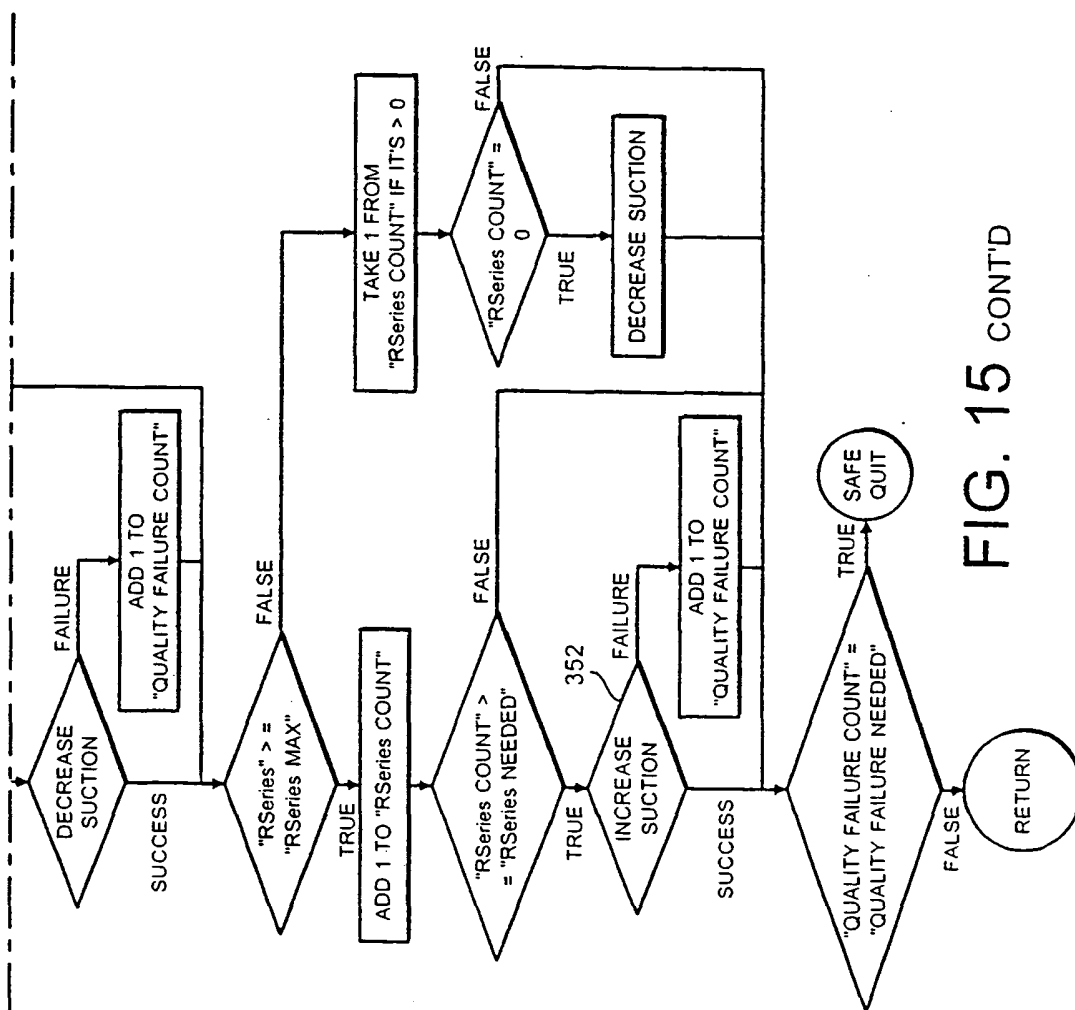


FIG. 15 CONT'D



FIG. 16

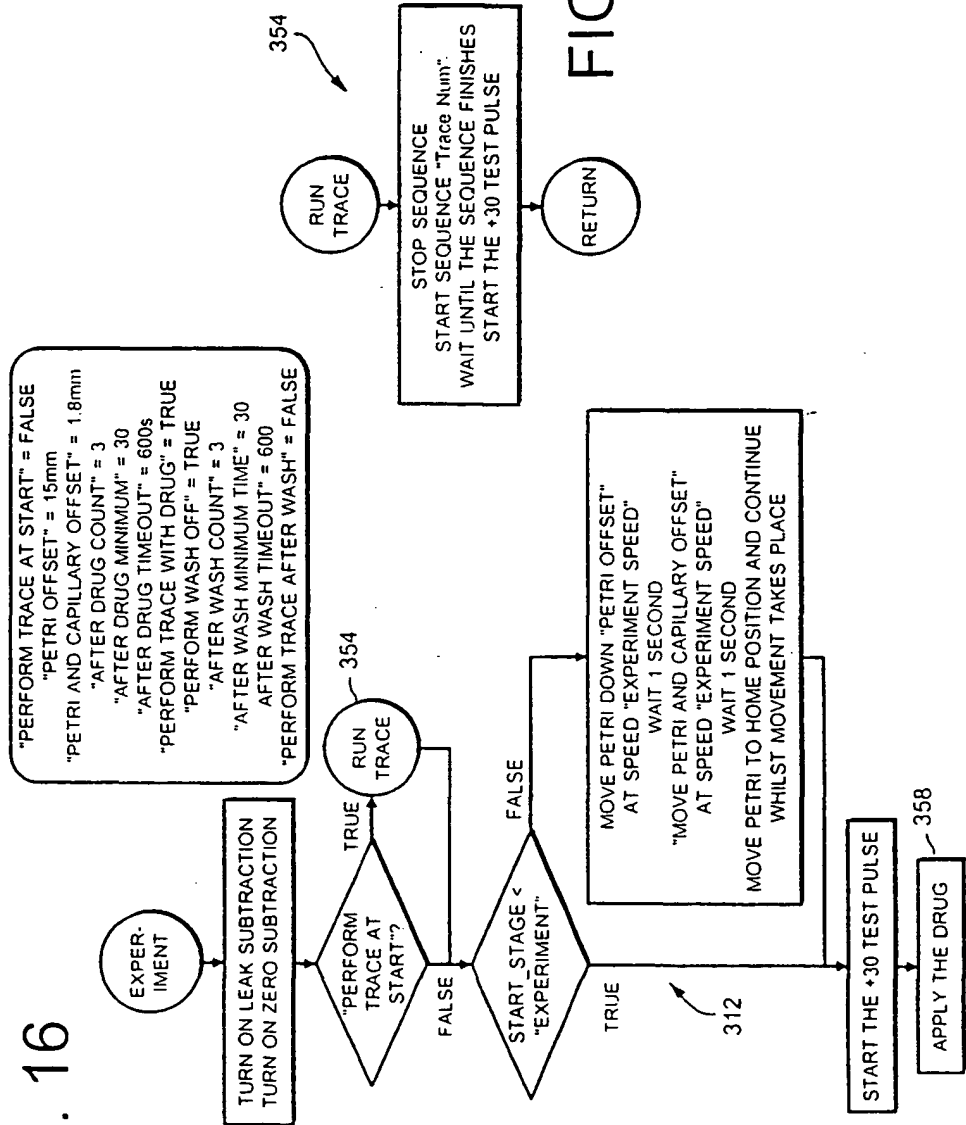

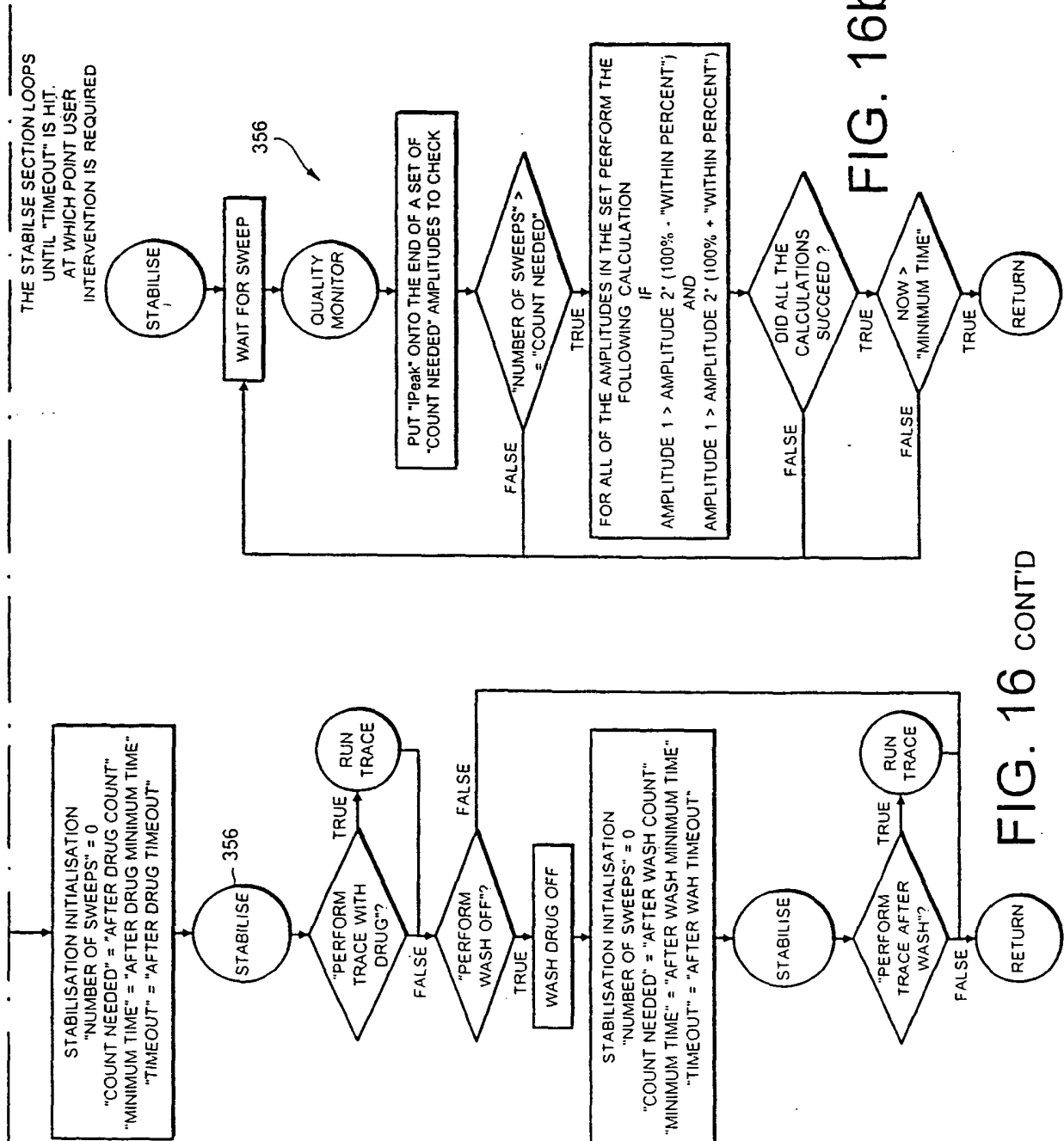

24 / 33


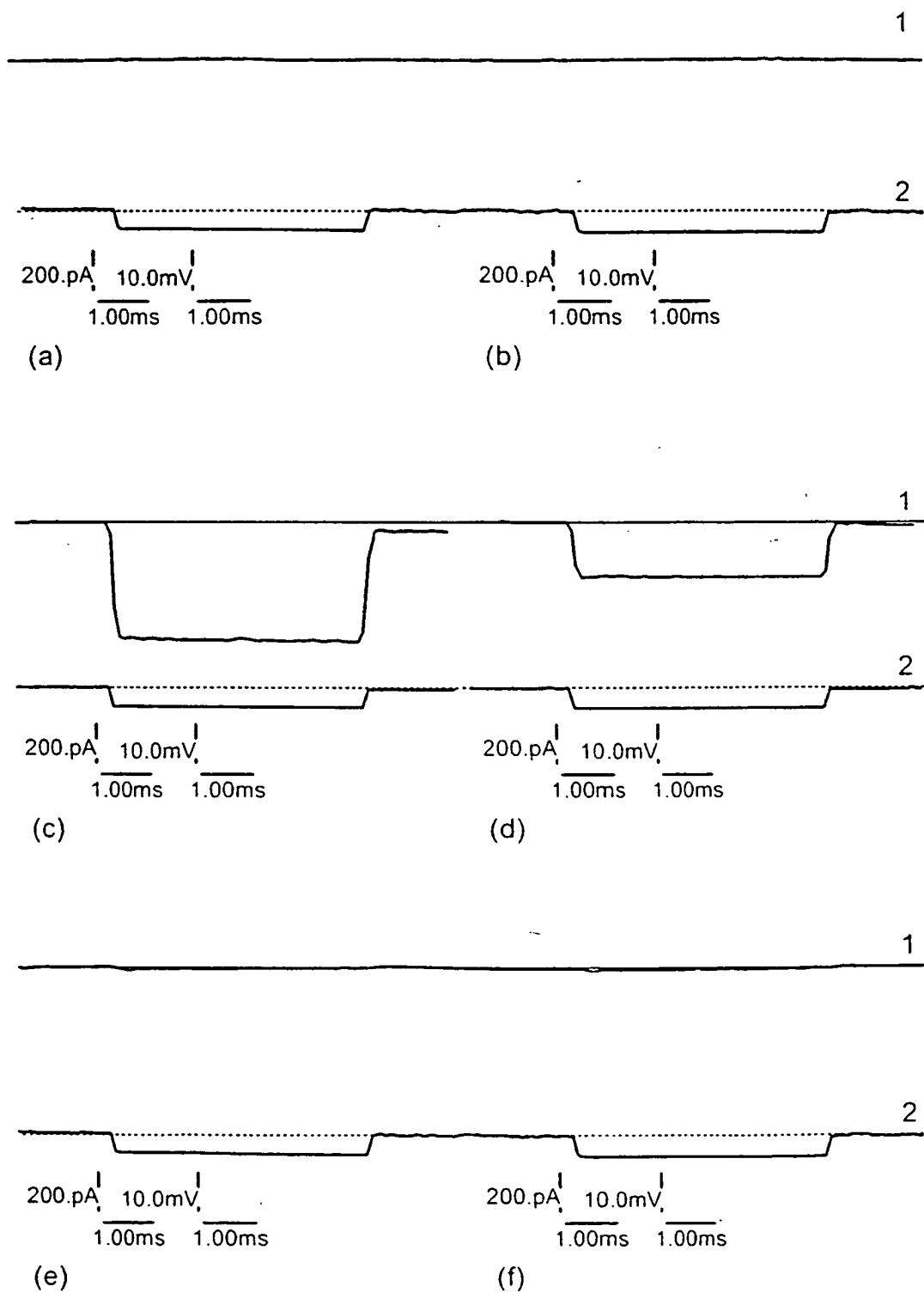


FIG. 17

25 / 33

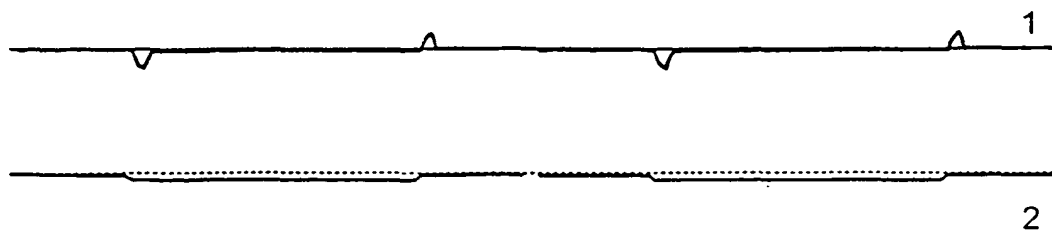


200.pA, 50.0mV,  
1.00ms 1.00ms

(a)

200.pA, 50.0mV,  
1.00ms 1.00ms

(b)

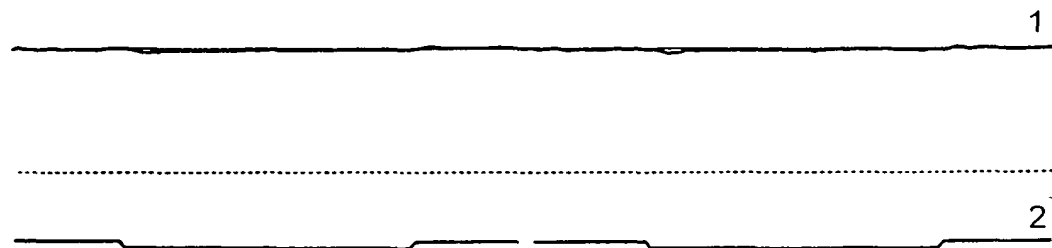

100.pA, 50.0mV,  
1.00ms 1.00ms

(c)

100.pA, 50.0mV,  
1.00ms 1.00ms

(d)

FIG. 18

26 / 33


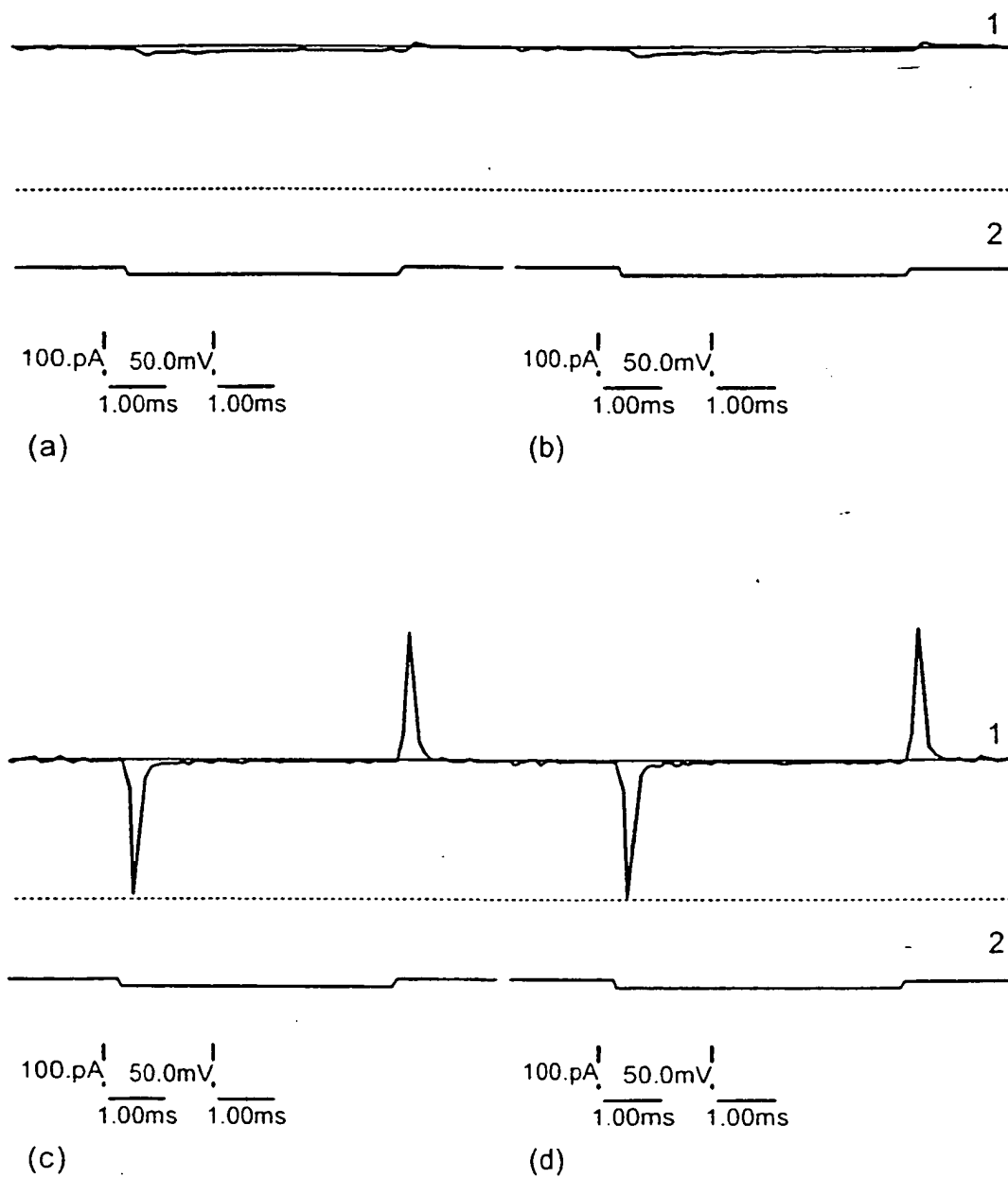



FIG. 19

27 / 33

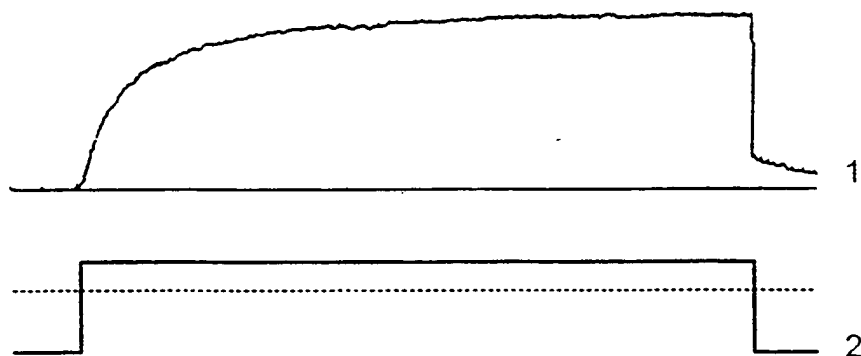

500.pA | 50.0mV |  
20.0ms 20.0ms

(a)

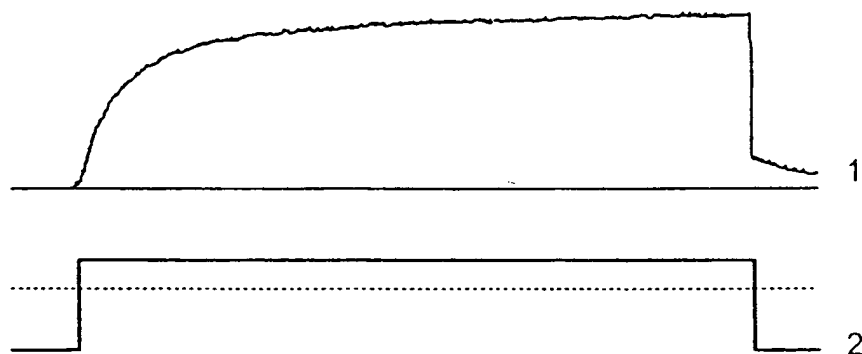

500.pA | 50.0mV |  
20.0ms 20.0ms

(b)


FIG. 20

28 / 33




200.pA |  
20.0ms

(a)

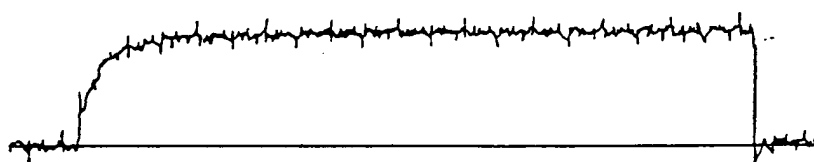

200.pA |  
20.0ms

(b)

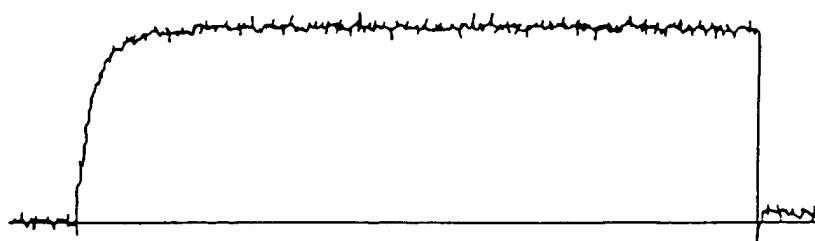

200.pA |  
20.0ms

(c)

FIG. 21

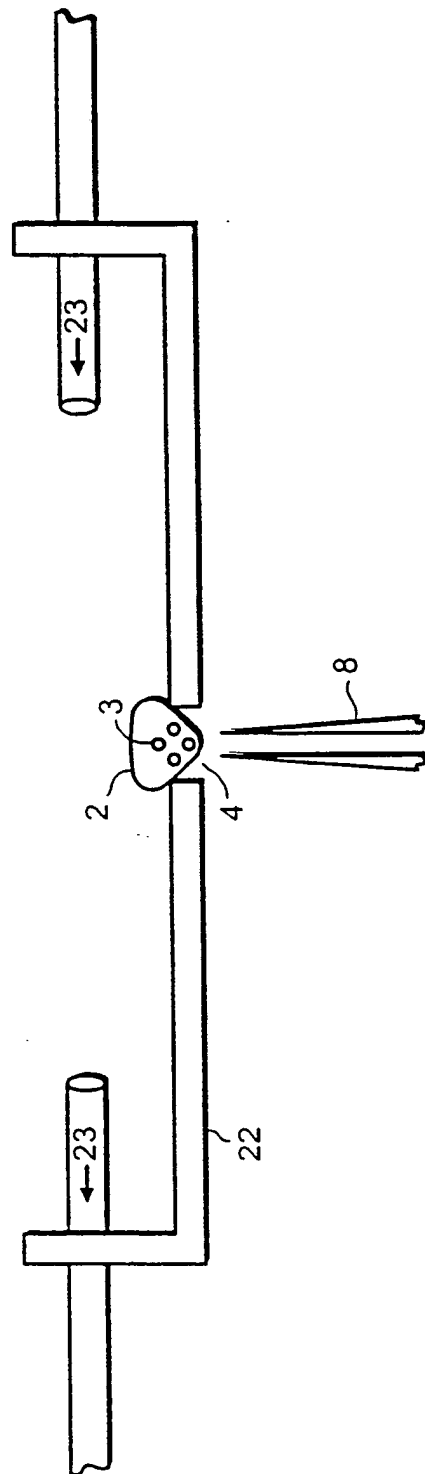
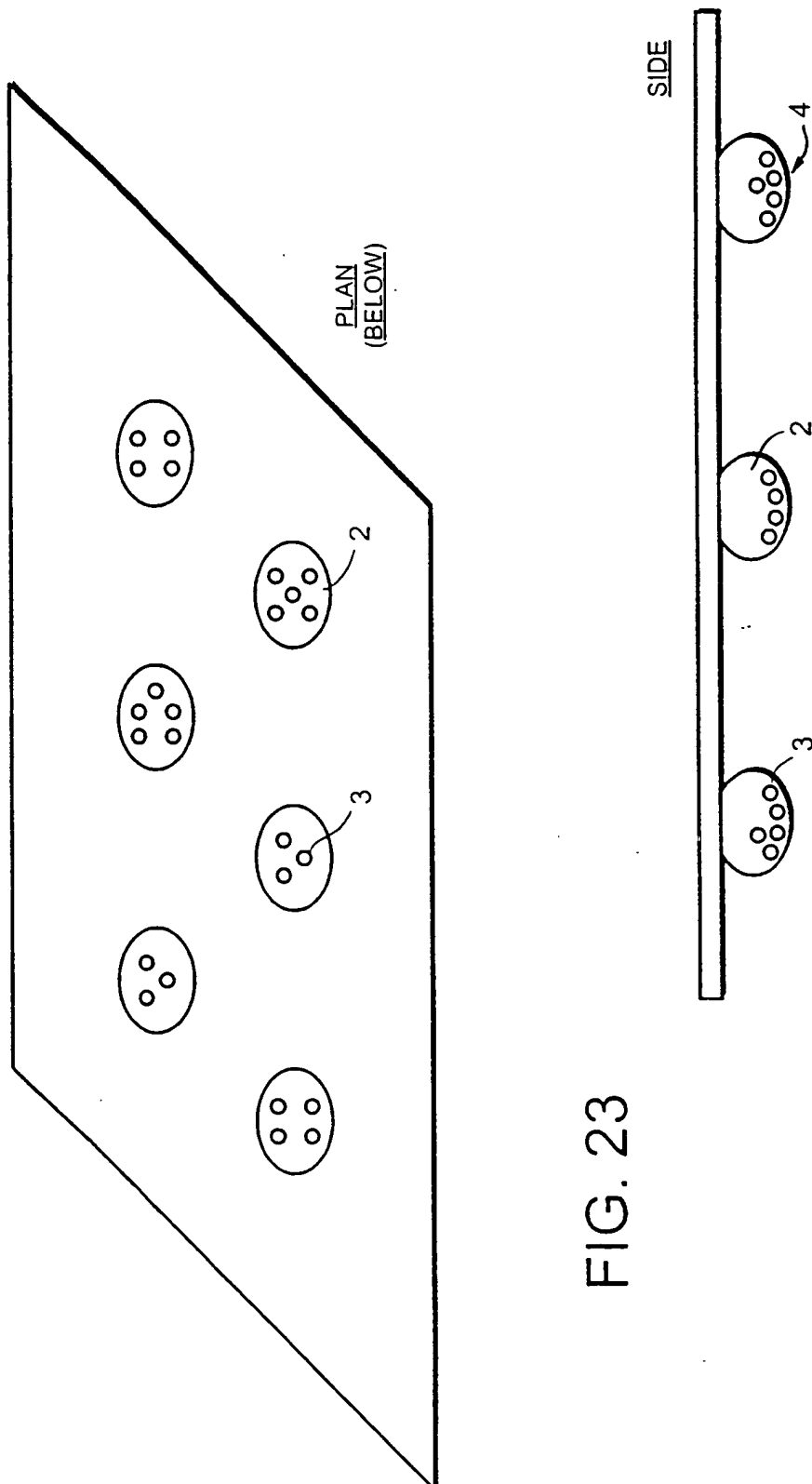


FIG. 22





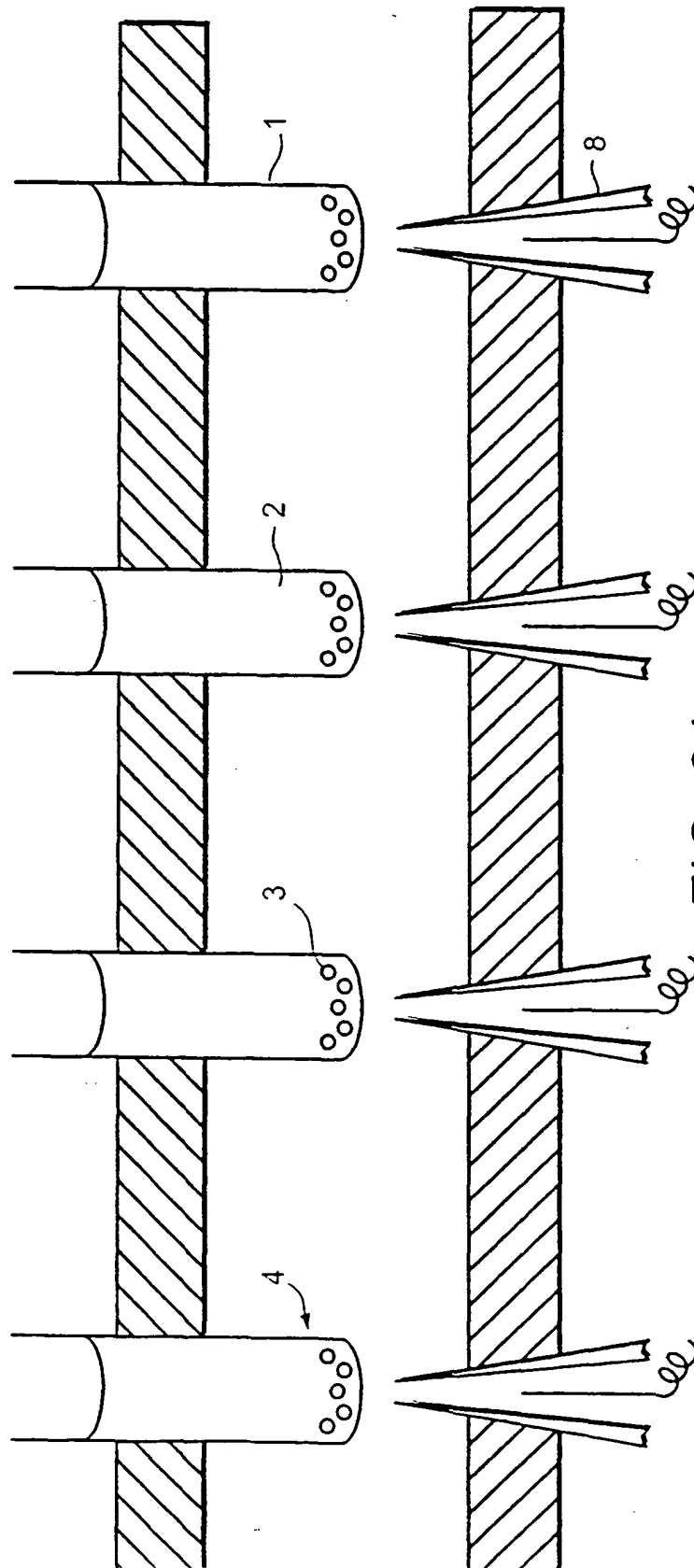


FIG. 24

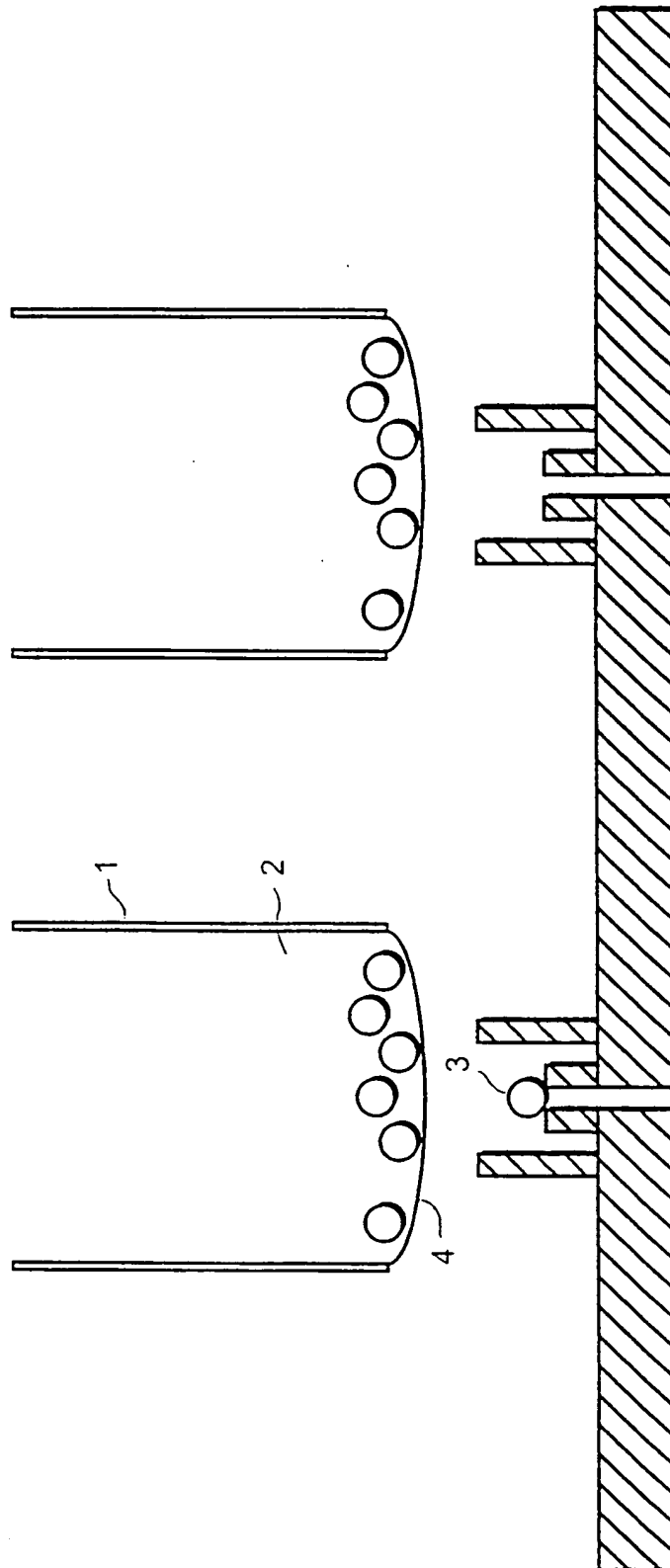


FIG. 25

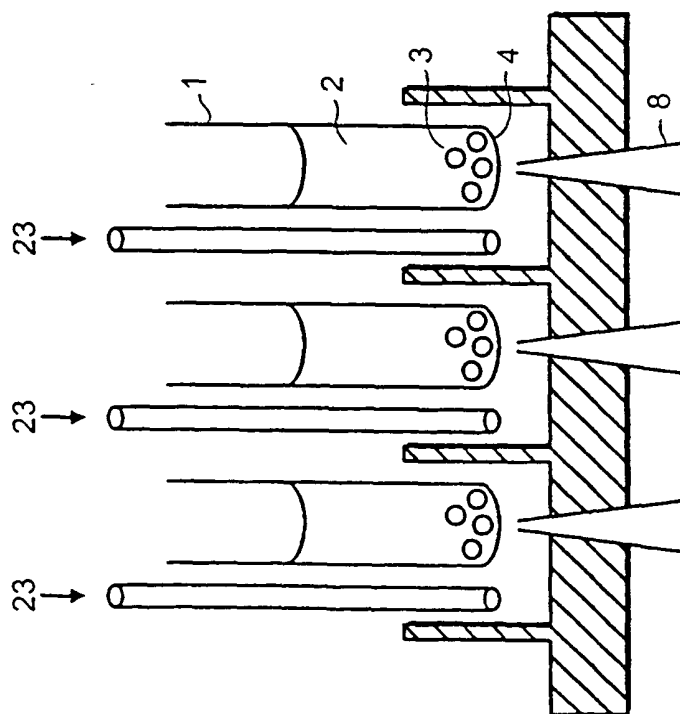


FIG. 26a

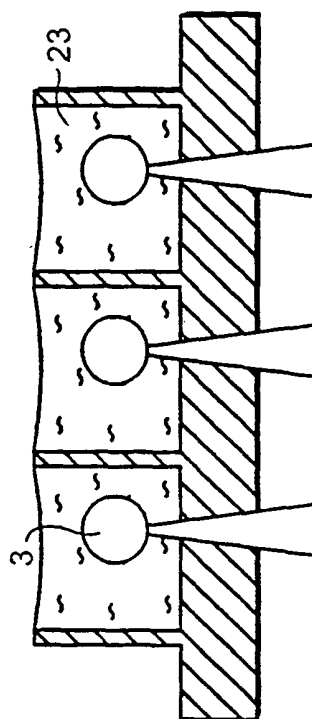


FIG. 26b

# INTERNATIONAL SEARCH REPORT

International Application No

PCT/GB 99/04073

## A. CLASSIFICATION OF SUBJECT MATTER

IPC 7 G01N33/487 G01N35/00 G01N35/10 //C12M1/34,C12N13/00

According to International Patent Classification (IPC) or to both national classification and IPC

## B. FIELDS SEARCHED

Minimum documentation searched (classification system followed by classification symbols)

IPC 7 G01N C12M

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the international search (name of data base and, where practical, search terms used)

## C. DOCUMENTS CONSIDERED TO BE RELEVANT

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
P, X	<p>DE 198 41 337 C (MICRONAS INTERMETALL GMBH) 23 September 1999 (1999-09-23)</p> <p>column 1, line 3 -column 1, line 43</p> <p>column 2, line 13 -column 2, line 31</p> <p>column 3, line 27 -column 3, line 44</p> <p>column 4, line 33 -column 4, line 62</p> <p>column 5, line 35 -column 5, line 41</p> <p>column 7, line 39 -column 7, line 51</p> <p>column 8, line 16 -column 9, line 48</p> <p>column 10, line 5 -column 10, line 28</p> <p>column 11, line 4 -column 11, line 23</p> <p>column 11, line 67 -column 12, line 35</p> <p>column 12, line 57 -column 13, line 22</p> <p>figures 1-9</p> <p style="text-align: center;">---</p> <p style="text-align: center;">-/--</p>	1-3,5-7, 9,12-14

☒ Further documents are listed in the continuation of box C.

☒ Patent family members are listed in annex.

### \* Special categories of cited documents :

"A" document defining the general state of the art which is not considered to be of particular relevance

"E" earlier document but published on or after the international filing date

"L" document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)

"O" document referring to an oral disclosure, use, exhibition or other means

"P" document published prior to the international filing date but later than the priority date claimed

"T" later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention

"X" document of particular relevance: the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone

"Y" document of particular relevance: the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art.

"&" document member of the same patent family

Date of the actual completion of the international search

28 February 2000

Date of mailing of the international search report

03/03/2000

Name and mailing address of the ISA

European Patent Office, P.B. 5818 Patentlaan 2  
NL - 2280 HV Rijswijk  
Tel. (+31-70) 340-2040, Tx. 31 651 epo nl,  
Fax: (+31-70) 340-2046

Authorized officer

Koch. A

# INTERNATIONAL SEARCH REPORT

Int. onal Application No

PCT/GB 99/04073

## C.(Continuation) DOCUMENTS CONSIDERED TO BE RELEVANT

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
X	<p>DE 197 12 309 A (NMI NATURWISSENSCHAFTLICHES UN) 20 May 1998 (1998-05-20) column 1, line 18 -column 1, line 43 column 1, line 63 -column 2, line 7 column 6, line 50 -column 7, line 49 column 8, line 8 -column 8, line 53 figures 1-7 column 3, line 46 -column 4, line 20 figures 1-3</p>	1,3,4
P,A	<p>DE 197 44 649 A (FRAUNHOFER GES FORSCHUNG) 15 April 1999 (1999-04-15) column 1, line 18 -column 1, line 43 column 2, line 17 -column 3, line 29 column 3, line 46 -column 4, line 20 figures 1-4</p>	1,9,11, 13,15-18
P,A	<p>WO 98 55870 A (A &amp; SCIENCE INVEST AB ;ORWAR OWE (SE); JARDEMARK KENT (SE)) 10 December 1998 (1998-12-10) page 1, line 2 -page 1, line 25 page 6, line 30 -page 7, line 5 page 20, line 22 -page 22, line 8 page 22, line 31 -page 23, line 24 page 23, line 36 -page 24, paragraph 1 figures 2-5</p>	1-3,12

# INTERNATIONAL SEARCH REPORT

Information on patent family members

International Application No

PC1/GB 99/04073

Patent document cited in search report	Publication date	Patent family member(s)	Publication date
DE 19841337 C	23-09-1999	EP 0962524 A JP 11346764 A DE 19827957 A EP 0960933 A JP 11346794 A	08-12-1999 21-12-1999 09-12-1999 01-12-1999 21-12-1999
DE 19712309 A	20-05-1998	WO 9822819 A EP 0938674 A	28-05-1998 01-09-1999
DE 19744649 A	15-04-1999	WO 9919729 A	22-04-1999
WO 9855870 A	10-12-1998	AU 8046898 A	21-12-1998

